

Cole-Parmer® Touch Screen

USER'S MANUAL



COLE-PARMER® TOUCH SCREEN SERIES



www.coleparmer.com

Table of Contents

SUBJECT	PAGE NO.
General Information.....	6
Safety Information	6
Warranty and Repair Information	7
Manual Description	7
Warranty	7
Repair Facilities and Parts	7
Serial Numbers	7
Calibrations	7
Specifications	8
Product Overview.....	9
Pump Models.....	10
Cole-Parmer® Touch Screen Series Button Reference.....	11
Selecting an Item from a Menu List	12
Messages	12
Entering Values	12
Operating Instructions.....	13
Power Connections	13
Syringe Loading.....	14
Configuring the Cole-Parmer® Touch Screen Series	15
Selection Mode	15
Syringe Selection.....	16
Flow Rate Selection	17
Target Volume/Time Selection	18
Previewing the Configuration	19

SUBJECT	PAGE NO.
Configuring the Pump's General Parameters	19
Setting the Time and Date	20
Setting the Audible Alarms.....	21
For Near End of Run Selection	22
Locking the Touchscreen	23
Setting the Force Limits	24
Setting the Display Brightness	25
Setting the Pump Address	26
Setting the RS 232 Baud Rate.....	26
Resume Run on Power Loss	27
Setting and Viewing Pump Information.....	27
Operating the Cole-Parmer® Touch Screen Series.....	28
Running the Pump	29
Custom Program Workflow	30
Step-by-Step Operating Example.....	31
Configuring Custom Programs	32
Setting a Time Delay for Program Operation.....	34
Creating a Custom New Program	35
Modifying an Existing/User-Defined Program	37
Changing the Order of Steps in a Program	37
Deleting a Step from a Program	37
Duplicate a Step in a Program	37
Pump Operation Profiles.....	38
Constant rate Profile	39
Ramp Profile	40
Stepped Profile	41
Pulse Profile	43

SUBJECT	PAGE NO.
Incorporating Advanced Options	44
Go To Step	44
Delay	46
Repeat Step.....	46
Link Program	47
Stop	47
I/O Event Triggers	48
I/O Trigger (OUT)	49
I/O Event (IN)	50
I/O On Event (IN).....	51
Managing Cole-Parmer® Touch Screen Series Programs	52
Saving a Copy of a Program	52
Renaming a Program	52
Appending a Program	53
Deleting a Program	53
Program Import and Export	54
Transfer a Program between Pumps	54
Exporting a Program to an External Computer	55
Importing a Program from an External Computer	55
External Pump Control	56
Setting up HyperTerminal	56
Pump Chain Commands.....	57
System Commands.....	58
Run Commands.....	59
Rate Commands	59
Volume Commands	60
Time Commands.....	60
Error Codes	61

SUBJECT	PAGE NO.
Appendices	63
Appendix A: Syringe Volume/Diameter Reference Table.....	63
Appendix B: Min/Max Flow Rates	64
Appendix C: Rack Accessories/Options	65
Mounting a 4x140 Multi-Rack	65
Mounting a 6/10 Multi-Rack.....	67
Syringe Loading – 4x140 & 6/10 Multi-Racks.....	69
Mounting a Micro Liter Multi-Rack	70
Syringe Loading – Micro Liter Multi-Rack	72
Syringe Loading – Push-Pull Mechanism	73
Appendix D: External Connections	75
User I/O Connector Specifications	75
Analog Input Option.....	76
Utilizing USB Interface.....	77
RS-232 Specifications	79
RS-485.....	80
Appendix E: Maintenance	81
Maintenance	81
Fuse Replacement	82
Battery Replacement	82
Upgrading Cole-Parmer® Touch Screen Series Software	83
Troubleshooting	84
Appendix F: 'Auto Fill' Valves.....	85
Setting Up the Valve	85
Appendix G: Ordering Information.....	86

General Information

SAFETY INFORMATION

Please read the following safety precautions to ensure proper use of your syringe pump. If the equipment is used in a manner not specified, the protection provided by the equipment may be impaired.

To Prevent Hazard or Injury:

USE PROPER LINE CORD

Use only the specified line cord for this product and make sure line cord is certified for country of use. The operating voltage range for the Cole-Parmer® Touch Screen Series is 100-240 vac, 50/60 Hz.

GROUND THE PRODUCT

This product is grounded through the grounding conductor of the power cord. To avoid electric shock, the grounding conductor must be connected to earth ground. Before making any connections to the input or output terminals of the product, ensure that the product is properly grounded.

MAKE PROPER CONNECTIONS

Make sure all connections are made properly and securely. Any signal wire connections to the unit must be no longer than 3 meters.

OBSERVE ALL TERMINAL RATINGS

Review the operating manual to learn the ratings on all connections.

USE PROPER FUSE

Use only specified fuses with product.

AVOID EXPOSED CIRCUITRY

Do not touch any electronic circuitry inside of the product.

AVOID PINCH HAZARD

A pinch hazard may exist between the pusher block and end blocks. Avoid placing fingers between these points while the pump is running.

DO NOT OPERATE WITH SUSPECTED FAILURES

If damage is suspected on or to the product do not operate the product. Contact qualified service personnel to perform inspection.

ORIENT THE EQUIPMENT PROPERLY

Do not orient the equipment so that it is difficult to operate the disconnection device.

PLACE PRODUCT IN PROPER ENVIRONMENT

Review the operating manual for guidelines for proper operating environments.

OBSERVE ALL WARNING LABELS ON PRODUCT

Read all labels on product to ensure proper usage.



CAUTION
Refer to Manual



Protective
Ground Terminal

WARRANTY AND REPAIR INFORMATION



CAUTION: REFER TO SAFETY INFORMATION AND SETTING UP THE Cole-Parmer® Touch Screen Series BEFORE PLUGGING IN THE PUMP.

Manual Description

This manual is designed to provide all operational and program information required to operate and maintain the Cole-Parmer® Touch Screen Series syringe pumps. The functions and features are described in the Technical Specifications section.

Warranty

The Cole-Parmer® Instrument Company warrants this product to be free from significant deviations in material and workmanship for a period of one year from date of purchase. If repair or adjustment is necessary and has not been the result of abuse or misuse within the one year period, please return-freight prepaid-and correction will be made with-out charge. Cole-Parmer® Instrument Company alone will determine if the product problem is due to deviations or customer misuse.

Out of warranty products will be repaired on a charge basis.

Return of Items

Authorization must be obtained from our Customer Satisfaction Department before returning items for any reason. When applying for authorization, please include data regarding the reason the items are to be returned.

For your protection, item must be carefully packed to prevent damage in shipment and insured against possible damage or loss. Cole-Parmer® Instrument Company will not be responsible for damage resulting from careless or insufficient packing. A restocking charge will be made on all unauthorized returns.

NOTE: the Cole-Parmer® Instrument Company reserves the right to make improvements in design, construction and appearance of our products without notice.

Repair Facilities and Parts

Cole-Parmer® stocks replacement and repair parts. When ordering, please describe parts as completely as possible, preferably using a part number obtained from our Customer Service department. If practical, enclose a sample part or sketch. We offer a complete reconditioning service.

Serial Numbers

All inquiries concerning our product should refer to the serial number of the unit, located on the rear panel.

Calibrations

All electrical apparatus are calibrated at rated voltage and frequency. While the flow and volume will stay calibrated, the peak pressure may vary. Cole-Parmer® recommends an annual calibration of the pump.

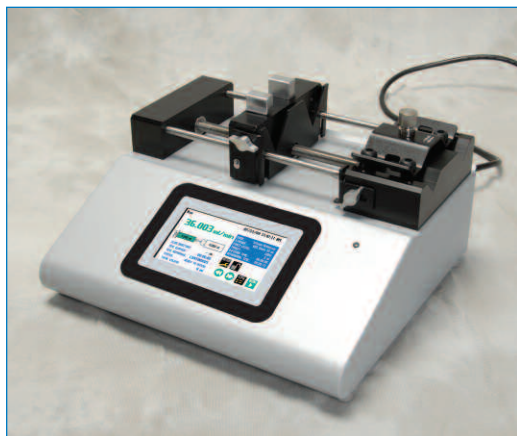


CAUTION: For research use only. Not for clinical use on patients.

SPECIFICATIONS

Accuracy	± 0.35%
Reproducibility	± 0.05%
Syringes (Min./Max.)	0.5 µl / 140 ml
Flow Rate:	
Minimum (0.5 µl syringe)	5 pl/min
Maximum (140 ml syringe)	215.803 ml/min
Display	4.3" WQVGA TFT Color Display with Touchpad
Non-Volatile Memory	Stores all settings
Connectors:	
RS-232	9 pin D-Sub Connector
RS-485	IEEE-1394, 6 pos
USB	Type B
I/O & TTL	15 pin D-Sub Connector
Linear Force (Max)	34 kg (75 lbs) @ 100% Force Selection
Drive Motor	1.8° Stepper Motor
Motor Drive Control	Microprocessor with 1/16 microstepping
Number of Microsteps per one rev. of Lead Screw	6,400
Step Rate:	
Minimum	27.5 sec/µstep
Maximum	26 µsec/µstep
Pusher Travel Rate:	
Minimum	0.36 µm/min
Maximum	190.80 mm/min
Power	100-240 VAC:50/60 Hz, 50 W, 0.5 A fuse
Dimensions	8.89 x 25.4 x 27.94 cm (3.5 x 10 x 11 in)
Weight	4.9 Kg (10.75 lb)
Atmospheric Specifications:	
Operating Temperature	4°C to 40°C (40°F to 104°F)
Storage Temperature	-10°C to 70°C (14°F to 158°F)
Humidity	20% to 80% RH, non condensing
Mode of Operation	Continuous
Classification	Class I
Pollution Degree	1
Installation Category	II
Supplier Name	Cole-Parmer®
Supplier Address	625 East Bunker Court , Verson Hills, IL 60061
Supplier Phone Number	800.323.4340
Regulatory Certifications	CE, UL, CSA, CB Scheme, EU RoHS

Product Overview



The Cole-Parmer® Touch Screen Series is the next generation of syringe pumps. The series offers unparalleled ease of use through the high resolution color touch screen user interface. The full touch screen interface enables the user to quickly create configurations and recall them for easy use. The TFT color display presents all the pump operating parameters on an easy to view run screen. The touch pad "lock" feature prevents accidental configuration changes ensuring the programs integrity is retained. The entire Cole-Parmer® Touch Screen Series features a battery backed real time clock. Programs can also be set to be started on a time delay controlled by a real time clock for unattended operation.

An led light on the front panel makes it easy to see if the pump is running. Once the pump starts to run, the display shows the pumps is running and the LED turns on. This will make it easy to see the pump's status.

The pumps are versatile and can be interconnected through the RS485 interface. Cole-Parmer® Touch Screen Series models can be mixed and matched in the daisy chain offering maximum flexibility. Up to 99 pumps can be linked together through the RS485 interface. This interface is easy to use provided each pump is assigned its own unique

pump address.



For limited laboratory space the Cole-Parmer® Touch Screen Series can be placed on its side to reduce the footprint by 55%.

The display also tilts with the change to allow the user operate the pump vertically.

Notifications and error messages are displayed for the user to acknowledge, eliminating any guesswork about problems.

Flow performance is optimized with a small step angle stepping motor that drives a precision lead screw and pusher block. The syringe mechanism is easy to use and securely holds the syringes for smooth flow performance.

Advanced microstepping techniques are employed to further reduce the step angle to eliminate flow pulsation. The accuracy is $\pm 0.35\%$ and 0.05% reproducibility. A wide dynamic flow range from 5 pl/min to 215.803 ml/min can be programmed into the pump. Flow rates are selectable with user selected engineering units from ml, μ l, nl pl, and hours, minutes and seconds.

Up to 34 kg (75 lbs) of linear force is available. This force is user adjustable to ensure the right force is applied for the various syringe sizes.

The superior design of the full metal chassis provides noise isolation and anti-vibration features for increased reliability. The spill dam is designed to prevent fluids from ingressing into the sealed display. All syringe racks are hardened rolled steel and will not deform with pressure.

Easy external connections to a computer or other control devices are made through USB interface or RS232 (9 pin Dsub). Simple ASCII commands make communication with the pump easy. For direct control of the pump, the user can use the I/O interface (15pin Dsub). Two trigger outputs and a trigger input for external events such as a process parameter is met. The footswitch input will allow the control of the pump through an external device. The unit also has an output for run indication allowing connection to a remote light.

All the models accommodate dual syringes except for the push/pull pumps which have two syringes on each side for simultaneous infusion and withdrawal. Complete flexibility is offered with the user able to use ANY manufacturers syringes with sizes ranging from 0.5 ul to 140 ml. Modular configuration is now done through the field installable syringe mechanisms. These optional accessories provide the necessary brackets to hold up to four 140 ml syringes or six 10 ul syringes. Also, the user can use the microliter rack which has been specifically designed for smaller syringes.

The programmable models offer maximum flexibility and capability for configuring and running different programs. Up to 40 programs of 20 steps each can be configured and stored in the unit, and recalled quickly with the touch of a button. The programs are identified with a 15 character alphanumeric name for easy identification and can be stored on the computer and downloaded at future dates. A real time clock is available to start and stop programs based on time. Extensive programming options allow the user to duplicate configurations, append one configuration to another, rename programs, delete programs, export programs to a PC or export a program to another pump. This feature can ensure consistency and accuracy of multiple pumps.

Preprogrammed profiles are selectable from Constant Rate, Ramp, Stepped and Pulse. These standard profiles make configuration easy to set up. If more complexity is needed the user can select from advanced programming options including:

- Go To statements
- Time delays
- Repeating steps
- Linking different programs
- Stopping the pump
- Triggering the pumps TTL output
- Accepting an event input in such as a user touch, motor stall
- Events can trigger the pump to withdraw or infuse

The pump Information/setup screen shows the pump's software version, calibration and lubrication intervals. The pump's alarm configurability includes alarms for near end of run (user selectable), completion of run, power up, keypad clicks, stall detection, calibration reminder.

Two options are available for the Cole-Parmer® Touch Screen Series. The analog input option which allows the analog control of the motor speed. By applying a 10 VDC max to the circuit the motor speed can be varied. The second option is for an internal fan. These will be factory installed.

Pump Models

The Cole-Parmer® Touch Screen Series is available in 5 different configurations.

Cole-Parmer® Touch Screen Infuse Only – Infuse only syringe pump. Accommodates 2 syringes. User definable flow rates with selectable target volumes or time values to control the total infusion volume.

Cole-Parmer® Touch Screen Infuse/Withdraw – Infuse/withdraw syringe pump. Accommodates 2 syringes. This unit support infuse only, withdraw only, infuse/withdraw, withdraw/infuse and continuous mode. User defined flow rates with selectable target volumes. One Program can be stored.

Cole-Parmer® Touch Screen Infuse/Withdraw Programmable – This is a Programmable Infuse/Withdraw with user defined configurations/programs of up to 800 steps. Up to 40 programs of 20 steps each can be stored in memory.

Cole-Parmer® Touch Screen Push/Pull – Push/Pull syringe pump. Accommodates 2 syringes for infusion and 2 syringes for withdrawal. This model supports infusion and withdrawal simultaneously at user defined flow rates and with selectable target volumes to control the total volume pumped. It also supports infuse only, withdraw only, infuse/withdraw, withdraw/infuse and continuous mode. User defined flow rates with selectable target volumes.

Cole-Parmer® Touch Screen Push/Pull Programmable – This is a Programmable Push/Pull pump with user defined program configurations of up to 800 steps. Multiple programs can be stored in memory.

Modular syringe racks can be purchased to create a multichannel syringe pump .

- 10 channel syringe rack
- 4 channel syringe rack
- µl syringe rack

COLE-PARMER® TOUCH SCREEN SERIES TOUCHSCREEN BUTTON REFERENCE

The following buttons are part of the Cole-Parmer® Touch Screen Series software:



Settings: Access the settings menu used to adjust the following parameters: Date and Time, Audible Alarms, Pusher Block Force, Display Brightness, RS-232 Communications, etc.



Fast Forward: Advances the Pusher Block in the forward (infuse) direction (press & hold).



Fast Reverse: Advances the Pusher Block in the reverse (withdraw) direction (press & hold).



Accept: Used to accept the parameters/settings on a screen and return to previous screen in the menu.



Cancel: Used to cancel any changes/entries on a screen and return to the previous screen in the menu.



Set Delay Parameters: Used to set the date and time values for delayed operation. This permits the pump to start automatically at a preset date and time.



Set Delay Parameters: Indicates the Delay Timer has been set. The pump will start automatically on the day/time selected.



Run Preview: Access the Run Preview screen that shows a summary operating parameters.



Run/Resume: Used to Run or Resume a Program.



Re-Run: Clears counters/timers (except the Total Infused Volume counter) and runs the selected program from the beginning.



Stop: Used to stop running a Program.



Reset: Resets all counters/timers.



Minimum Value Allowed: Used to select minimum value allowed.



Maximum Value Allowed: Used to select maximum value allowed.



Configure: Used to set up the program.



Back: Used to call the Mode and Programs menu.



Page Down: Used to page down in a display list.



Page Up: Used to page up in a display list.



Line Up: Used to scroll up in a display list.



Line Down: Used to scroll down in a display list.



Insert Step: Used to insert a new Step in a Program while editing/creating the Program. The new Step is inserted after the Step that is currently selected.



Delete Step: Used to delete the selected Step in a Program while editing/creating the Program. To confirm deletion, you must select the Delete Step button a second time.



Move Step Up: Moves the currently selected Step in a Program up one position in the Program Order (e.g. Step 3 becomes Step 2).



Move Step Down: Moves the currently selected Step in a Program down one position in the Program Order (e.g. Step 2 becomes Step 3).



Copy Step: Creates a duplicate (copy) of the currently selected Step in a Program and inserts it as the last Step in the Program being created/edited.



Touchscreen Unlocked: Indicates that the Touchpad is currently unlocked when the Touchpad Lock Enabled parameter is selected in the Settings menu. Pressing the button locks the Touchpad.



Touchscreen Locked: Indicates that the Touchpad is currently locked when the Touchpad Lock Enabled parameter is selected in the Settings menu. Press the button twice to unlock the Touchpad.



File Options: Access the File Options menu that is used to manage the list of Programs and to export Program information to an externally linked pump or to a connected computer.

Selecting An Item From a Menu List

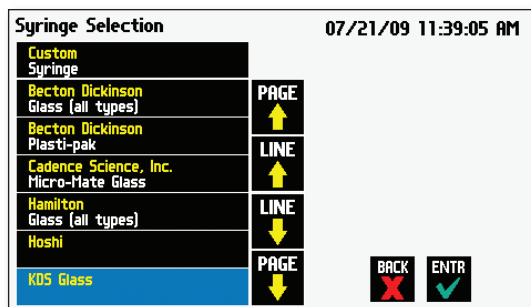
Selecting items is easy with the Cole-Parmer® Touch Screen Series.

Touch the item to select/highlight the item, then touch again to select it or press **ENTR**



Canceling a Selection

Canceling selections is easy.





Messages

There are different types of messages. Error messages are shown in yellow.

- Critical
- Run-mode
- Program
- Set-up

Notices, alarms, confirming, and verified messages are shown in green.

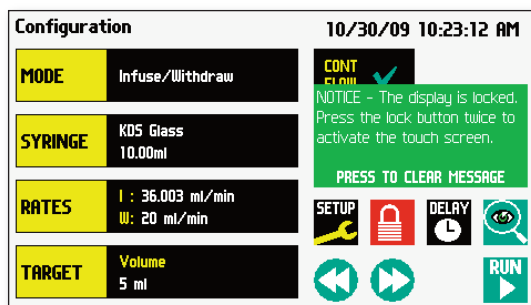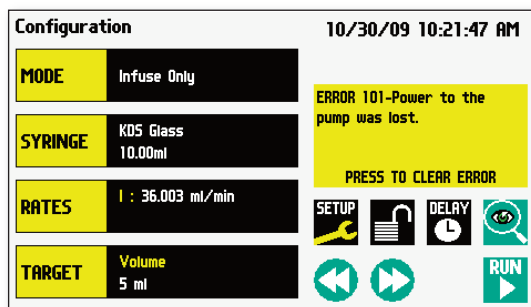- Alarms
- Confirm
- Verified



Entering Values

Entering new values is easy.

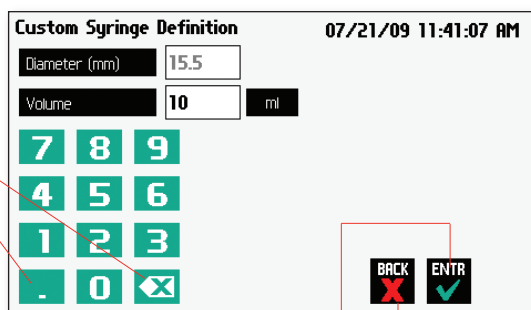
Use the keypad to enter the number

Backspace Key


Decimal point

Accept the value and go back up to the previous screen.

Cancel and go back to the previous screen.



OPERATING INSTRUCTIONS

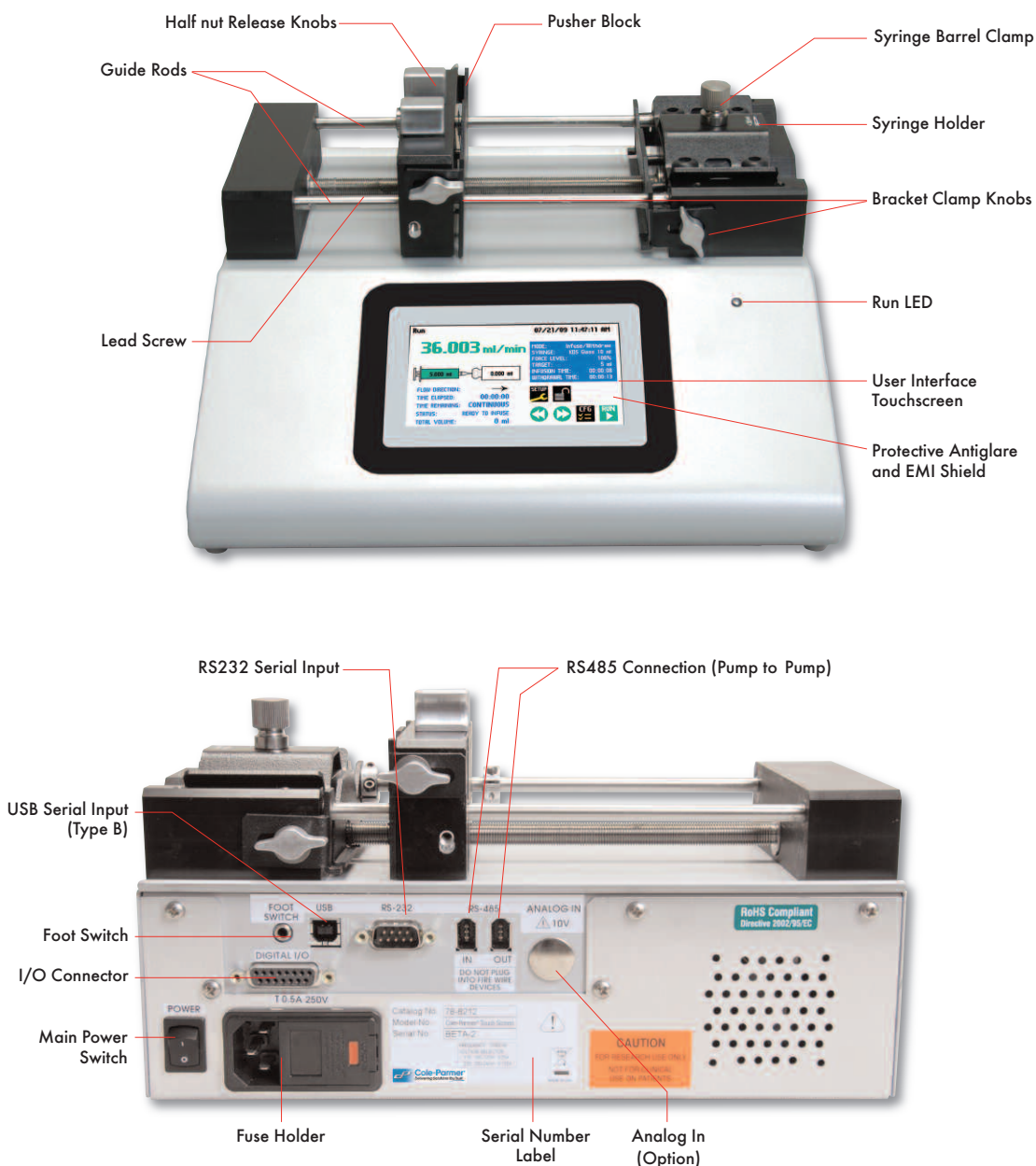

POWER CONNECTIONS

1. Plug the power cord in. The Cole-Parmer® Touch Screen Series has a universal power input connector on the rear of the unit. The operating range for the Cole-Parmer® Touch Screen Series series is 100 – 240 VAC, 50/60 Hz.
2. Turn on the main power switch.

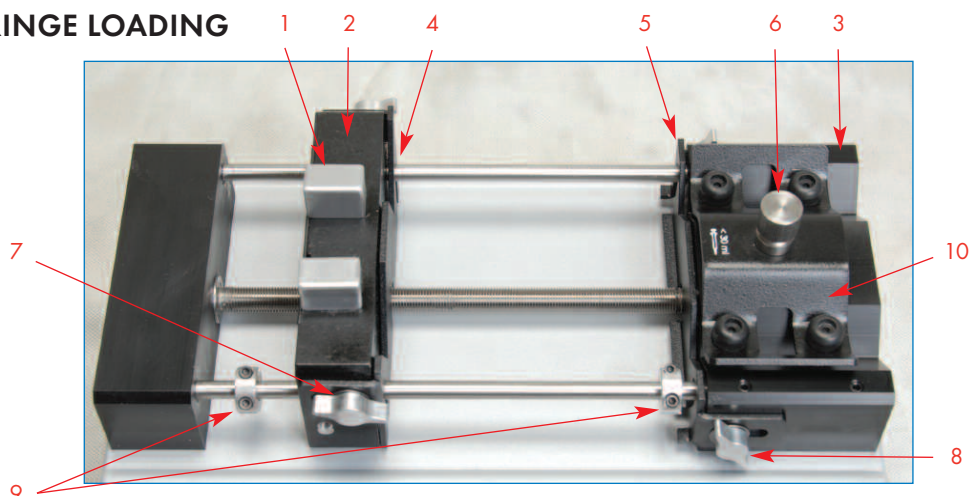


CAUTION: Do not connect to firewire ports on a computer. Damage may occur to the pump or computer.



NOTE: Remove tab under battery cover to activate battery back-up of date and time in the pump.

SYRINGE LOADING



The Standard Cole-Parmer® Touch Screen Series can accommodate up to two syringes. Any type of syringe, glass, plastic or stainless steel can be used. The pump can accommodate from 0.5 ul to 140 ml syringes. The syringe barrel clamp should be flipped for syringes smaller than 30 ml. For syringes larger the 30 ml the syringe clamp should be used as shown.

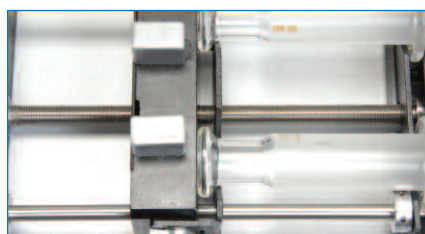
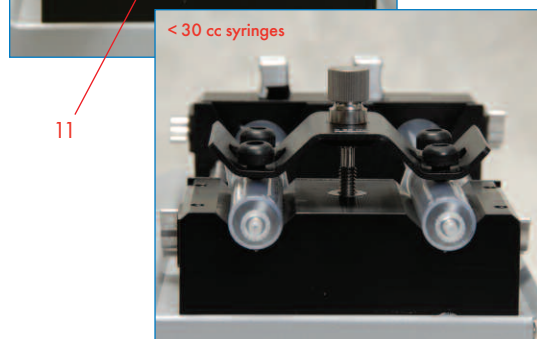
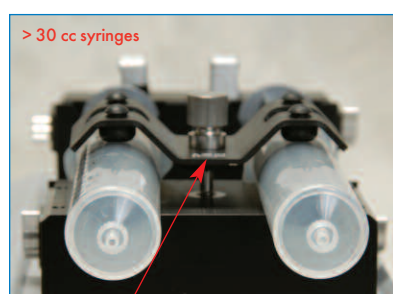
1. Squeeze the release tabs (1) and move the pusher block (2) to the approximate length of the syringe.
2. Unscrew the bracket clamping knobs (7,8) to loosen the retaining brackets (4,5) on the pusher block (2) and the Syringe holder block (3).
3. Unscrew the syringe barrel clamp by unscrewing the syringe barrel clamp knob (6).
4. Place the syringe on the syringe holder block (3). Ensure the barrel flange is with in the syringe holder's block's retaining bracket (5) and the plunger flange is

in the pusher' blocks retaining bracket.

5. Screw down the syringe clamp knob (6) until the syringe barrel clamp (10) is tight against the syringe barrel. Do not overtighten, especially when using glass syringes.
6. Secure the syringe flanges by pushing the retaining brackets (4,5) tightly against the syringe flanges using the bracket clamping knobs (7,8). The syringe plunger should be secured by the pusher block retaining bracket and the flange of the syringe should be secured by the syringe holder block retaining bracket.
7. Tighten the mechanical stop collars (9) on the pump to the mechanical stop for refill and infusion.



NOTE: The syringe barrel clamp (10) may need to be flipped depending on your syringe size. To do so, completely unscrew and remove the syringe barrel clamp knob from the syringe block. Then remove lock nut (11). Flip the bracket and reassemble.



CAUTION: Be sure mechanical stop collars are positioned properly to prevent the plunger from bottoming out on glass syringes.

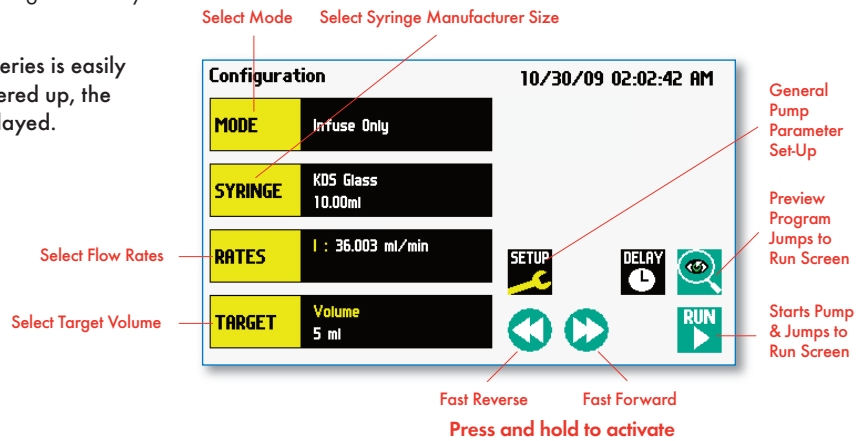


NOTE: For infuse/withdraw models make sure the plunger is secured in the retaining bracket and the barrel of the syringe is secured in the retaining bracket.

Configuring the Cole-Parmer® Touch Screen Series

The Cole-Parmer® Touch Screen Series program is easy to set up.


The Cole-Parmer® Touch Screen Series is easily configured. When the unit is powered up, the Configuration Screen will be displayed.



If you select infuse/withdraw or withdraw/infuse mode the configuration screen will have the continuous mode selection. (Press to toggle on/off)

Infuse/Withdraw models only. Continuous flow is an option on the display. In continuous mode, the pump will infuse then withdraw continuously or vice versa depending on which mode is selected. The pump will switch to withdrawal based on the target volume entered.

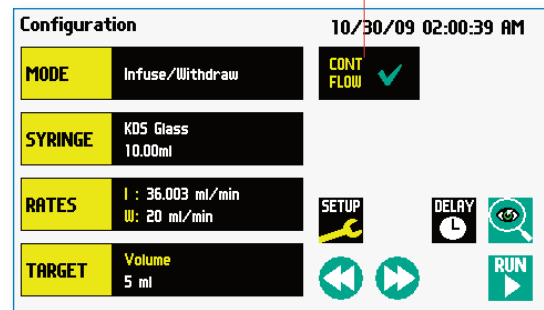The total infuse volume will be displayed on the run screen as well as total elapsed time.



Mode Selection

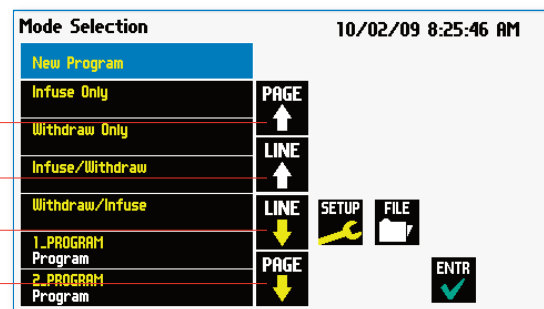
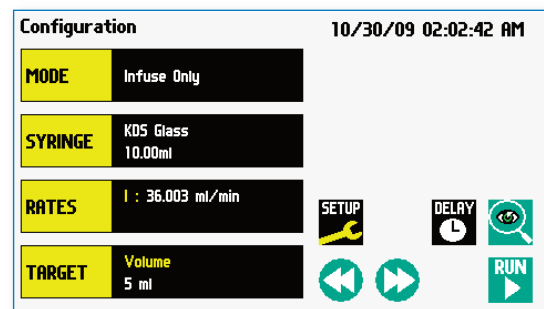
In the Configuration Display select mode (available modes depend on the pump model)

- Infuse
- Withdraw
- Infuse/Withdraw
- Infuse/Withdraw (continuous)
- Withdraw/Infuse
- Withdraw/Infuse (continuous)
- Custom program

The Mode Selection Screen will be displayed. The available modes will be displayed based on the model of syringe pump.

Press the Mode to select

- View programs up a page
- Views up one program
- Views down one program
- View the programs down a page



Syringe Selection

1. In the **Configuration Display** select the syringe button.

Configuration 10/30/09 02:05:11 AM

MODE	Infuse Only
SYRINGE	KDS Glass 10.00ml
RATES	I : 36.003 ml/min
TARGET	Volume 5 ml

SETUP DELAY RUN

2. Select the manufacturer type. Then the available syringe sizes will be displayed. In this example we have selected a KDS 10ml Glass Syringe.

Syringe Selection 07/21/09 11:39:05 AM

Custom Syringe	PAGE
Becton Dickinson Glass (all types)	LINE
Becton Dickinson Plasti-pak	LINE
Cadence Science, Inc. Micro-Mate Glass	LINE
Hamilton Glass (all types)	LINE
Hoshi	PAGE
KDS Glass	PAGE

BACK ENTR

Press

Syringe Size Selection 07/21/09 11:39:46 AM

KDS Glass

1 ml	4.8 mm	PAGE
2 ml	6.45 mm	LINE
5 ml	12.6 mm	LINE
10 ml	15.5 mm	LINE
20 ml	20.4 mm	PAGE
30 ml	22.9 mm	PAGE

BACK ENTR

Press

3. For a Custom Syringe use this selection if syringe manufacturer is not listed.

Press to enter a custom syringe

Syringe Selection 07/21/09 11:39:05 AM

Custom Syringe	PAGE
Becton Dickinson Glass (all types)	LINE
Becton Dickinson Plasti-pak	LINE
Cadence Science, Inc. Micro-Mate Glass	LINE
Hamilton Glass (all types)	LINE
Hoshi	PAGE
KDS Glass	PAGE

BACK ENTR

Custom Syringe Definition Display

Press to enter the total syringe volume

Press to enter the inside diameter of the syringe

Custom Syringe Definition 07/21/09 11:41:07 AM

Diameter (mm)	15.5
Volume	10 ml

7 8 9
4 5 6
1 2 3
. 0 X

BACK ENTR

Press to change the units (µl, ml)



CAUTION: BE SURE TO LOWER THE FORCE VALUE WHEN USING GLASS SYRINGES.

Flow Rate Selection

1. In the **Configuration Display** enter the flow rate(s). In this example we want to enter 36.003 ml/min.

Press

The Configuration Display shows the following settings:

- MODE:** Infuse Only
- SYRINGE:** KDS Glass 10.00ml
- RATES:** I : 36.003 ml/min
- TARGET:** Volume 5 ml

Navigation buttons include SETUP, DELAY, and RUN. Arrows indicate the flow rate is being entered.

2. Depending on the Mode, the appropriate set rate screen will be displayed. To change the units, press the units.

Enter the flow rate

Press to change the units

The Set Infuse Rate screen shows the following settings:

- Infuse Rate:** 36.003 ml/min
- MAX:** 36.003 ml/min
- MIN:** 69.341 nl/min

Navigation buttons include BACK and ENTR. Arrows indicate the unit is being changed.

Press to select the maximum volume for the syringe size selected

Press to select the minimum volume for the syringe size selected

3. In the **Units Selection Display** select the volume and the time base.

The Units Selection screen shows the following settings:

- Current selection:** ml/min
- ml** (selected)
- ul**
- nl**
- pl**
- Hr**
- Min**
- Sec**

Navigation buttons include BACK and ENTR.

Target Volume/Time Selection

1. In the **Configuration Display** enter the target volume or time. In this example we want a volume of 5 ml.



NOTE: Target time is only allowed in infusion only or withdrawal only modes.

2. Depending on the Mode either Target Volume (I/W or W/I) will appear or Target Volume or Target Time Selection (I only or W only).



NOTE: For time enter h:m:s. or press again to select secs. Time is accurate to +/- 0.01 secs. Remember to enter the :

3. To change the units, press the units.

Press to select target volume

Press to select target time

Press to select the maximum volume for the syringe size selected

Configuration 10/30/09 02:09:04 AM

MODE Infuse Only

SYRINGE KDS Glass 10.00ml

RATES I : 36.003 ml/min

TARGET Volume 5 ml

SETUP DELAY RUN

Press Here

Press to change the units

Set Target Volume or Time 07/21/09 12:01:51 PM

Volume 5 ml

Time 0 secs

7 8 9 MAX 10 ml

4 5 6

1 2 3

. 0 X

BACK X ENTR ✓

4. In the **Units Selection Display** select the volume units.

Units Selection 06/25/09 8:21:56 AM

Current selection: ml

ml

ul

nl

pl


BACK X ENTR ✓

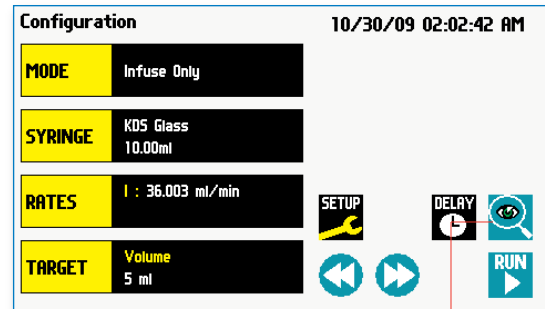
Previewing the Configuration

It is simple to view the Pump Configuration on the RUN display.

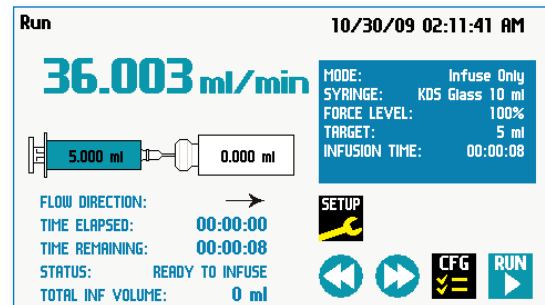
Press  to view Pump Configuration

After you have verified the configuration.

Press  to run the displayed configuration.




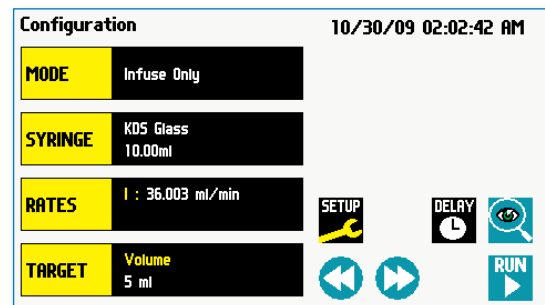
Press to preview configuration



Configuring the Pump's General Parameters

The pump has basic parameters that should be configured during initial start-up. Once these parameters are entered they probably won't be changed.

To enter the set-up screen press 



Setting the Time and Date

The date and time is displayed on each screen. A built in battery backup will maintain the date and time even if the pump loses power or is unplugged for a period of time.

In the Set-up screen, press **Date/Time**

Press



NOTE: The date AND time must be entered.



NOTE: Remember to activate the battery by removing the plastic tab on the bottom of the pump.

Set-up 08/24/09 9:17:18 AM

Date/Time	
Audible Alarms	Pump Info
Force 100%	Baud Rate 115200
Display Brightness 100%	Pump Address 0
Touchpad Lock Enable	Resume Run After Pwr Loss

BACK ENTR

Date/Time Display

Press **mm/dd/yy** to change the date format

- (month/day/year)
- (day/month/year)
- (year/month/day)

Enter the date.



NOTE: Remember to use the /



NOTE: The year must be entered as YY, i.e. 2009 is entered as 09

Press **AM/PM** to change the time format.

- AM
- PM
- 24 Hour

Enter the current time in hh:mm:ss (hours:minutes:seconds).



NOTE: Remember to use the :

Accept the changes made by pressing and go back to the Set-up screen.



Cancel the changes made by pressing and go back to the Set-up screen.



Set Date and Time 07/23/09 8:42:56 AM

Date: 07/23/09 mm/dd/yy

Time: 08:42:50 AM

7 8 9

4 5 6

1 2 3

/ 0 [X]

BACK ENTR

Press to change the date format

Press to change the time format

Setting the Audible Alarms

In the Set-up screen, press **Audible Alarms**.

Select the alarms, a check mark will indicate the alarms selected.

If these alarms are selected the pump will provide an audible alarm.

Accept the changes made by pressing **ENTR** and go back to the Set-up screen.

Cancel the changes made by pressing **BACK** and go back to the Set-up screen.

Set-up 08/24/09 9:17:18 AM

Date/Time	Upgrade Software
Audible Alarms	Pump Info
Force 100%	Baud Rate 115200
Display Brightness 100%	Pump Address 0
Touchpad Lock Enable	Resume Run After Pwr Loss

BACK ENTR

End of Run – Alarm when the program is finished

Power Up – Pump is turned on the pump

Keypad Clicks – When the keypad is touched the unit will provide audible feedback

Audible Alarms 08/24/09 9:18:41 AM

✓ End of Run	Power Up
Stall	Keypad Clicks
	✓ Calibration Reminder
Near End of Run	0%
Mute All	

BACK ENTR

Stall – Motor stalls

Near End of Run – When the volume dispensed is reached for the specified %

Mute All – No audible sounds

Calibration Reminder – Alarm when the unit has reached its calibration due date



NOTE: If power is removed from the pump without turning the pump switch off an audible alarm will sound. To stop the audible alarm toggle the power switch to off.

For Near End of Run Selection

Select the % volume dispensed. For example if you select 80%. The pump will sound an alarm (3 short beeps) when the 80% of the volume is dispensed.

Accept the changes made by pressing **ENTR** and go back to the Set-up screen.



Cancel the changes made by pressing **BACK** and go back to the Set-up screen.



Press to select % volume dispensed

For a Custom %

Use this selection if the % is not listed.

Press **Other Value** to enter a custom %

Accept the changes made by pressing **ENTR** and go back to the Set-up screen.



Cancel the changes made by pressing **BACK** and go back to the Set-up screen.



Press to enter a custom %

Enter the percentage.


Accept the changes made by pressing **ENTR** and go back to the Set-up screen.



Cancel the changes made by pressing **BACK** and go back to the Set-up screen.




Locking the Touchscreen


A touchscreen lock function can be enabled to prevent accidental operation or interruption of operation. The touch screen lock function makes all buttons on the Run display inactive except for the  button.



NOTE: This only enables the function, the user needs to activate it in the Run screen or Configuration screen.

In the **Set-up Display** select **Touchpad Lock Enable** function.

Accept the changes made by pressing  and go back to the Set-up screen.


Cancel the changes made by pressing  and go back to the Set-up screen.

Touchscreen Lock On

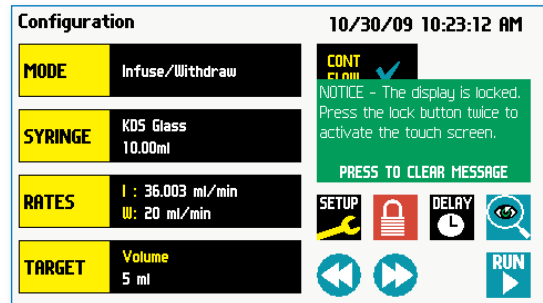
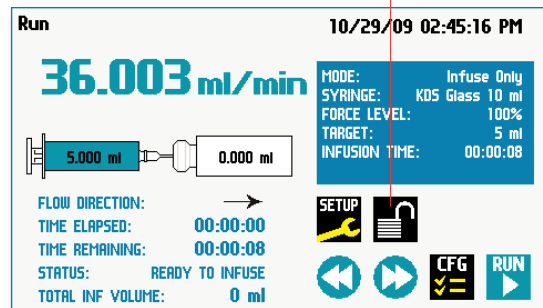
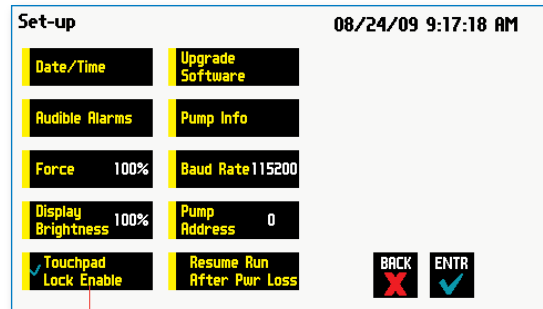
To turn the touchscreen lock on, make sure that the touchscreen lock is enabled. If it is, the button will be shown on the Run display.

Press the  it will turn .

Touchscreen Lock Off

In the Run display touch the  to turn the touchscreen lock off.

A Notice will ask you to touch the lock button again.



NOTE: The Stop button will always remain active allowing you to manually pause/stop pump operation even when the touchscreen is locked.

Setting the Force Limits

The maximum linear force of the syringe pump is 75 lbs or 35 kgs. This force is sufficient to damage, delicate, low volume syringes, as well as many standard syringes if pumping viscous fluids. The % of maximum force applied by the pusher block can be changed.

In the set up screen, press **Force**.



NOTE: A minimum value of approximately 30% is required to overcome friction forces in the mechanism. You may have to set the force level to a higher value depending on the backpressure/friction in your system.

Set the force level.

Press **Other Value** to enter a custom force %

Accept the changes made by pressing **ENTR** and go back to the Set Force screen.

Cancel the changes made by pressing **BACK** and go back to the Set Force screen.

For a Custom Force %

Enter the percentage for custom force %

Accept the changes made by pressing **ENTR** and go back to the Set Force screen.

Cancel the changes made by pressing **BACK** and go back to the Set Force screen.

Set-up 08/24/09 9:17:18 AM

Date/Time	Upgrade Software
Audible Alarms	Pump Info
Force 100%	Baud Rate 115200
Display Brightness 100%	Pump Address 0
Touchpad Lock Enable	Resume Run After Pwr Loss

BACK ENTR

Press Here

Set Force Level 07/21/09 11:27:59 AM

100%	50%
90%	40%
80%	30%
70%	20%
60%	Other Value

BACK ENTR

Press here to enter custom force %

Set Force Level 08/24/09 9:31:24 AM

Force 100

7	8	9
4	5	6
1	2	3
0	X	

BACK ENTR

Suggested Force Level Settings for Common Syringes

Syringe Type/Material	Capacity	Force Setting	Common Manufacturer(s)
Plastic Syringes	≤20ml	50%	BD plastic, Sherwood, Airtite, Terumo
	>20ml	100%	
Stainless Steel	All	100%	KDS, Cole-Parmer
Glass/Glass	≤20ml	20%	KDS Glass, Cadence (Popper)
	>20ml	30%	
Glass/Plastic	≤1000μl	30%	Hamilton, SGE
	≤5ml	50%	
	>5ml	100%	

Setting Display Brightness

The backlight can be adjusted.

In the Set-up screen, select **Display Brightness**

Press to select % Display Brightness

Select the % Display Brightness. For example if you select 80%. The brightness of the display will be 80% of the maximum.

Accept the changes made by pressing **ENTR** and go back to the Set-up screen.

Cancel the changes made by pressing **BACK** and go back to the Set-up screen.

For a Custom Value

Use this selection if the % is not listed.

Press **Other Value** to enter a custom %

Enter the percentage.

Accept the changes made by pressing **ENTR** and go back to the Set-up screen.

Cancel the changes made by pressing **BACK** and go back to the Set-up screen.

Set-up 08/24/09 9:17:18 AM

Date/Time	Upgrade Software
Audible Alarms	Pump Info
Force 100%	Baud Rate 115200
Display Brightness 100%	Pump Address 0
Touchpad Lock Enable	Resume Run After Pwr Loss

BACK ENTR

Press Here

Set Display Brightness 07/21/09 11:31:47 AM

100%	50%
90%	40%
80%	30%
70%	20%
60%	Other Value

BACK ENTR

Press here for custom display brightness level

Set Display Brightness 08/24/09 9:32:44 AM

Backlight 45

7	8	9
4	5	6
1	2	3
0	←	

BACK ENTR

Setting the Pump Address

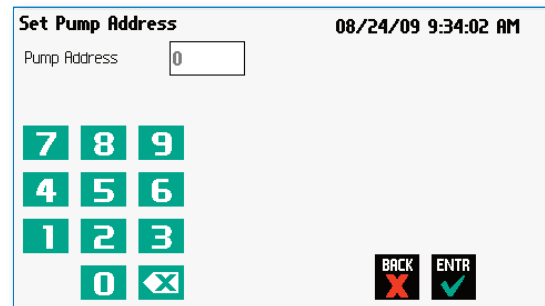
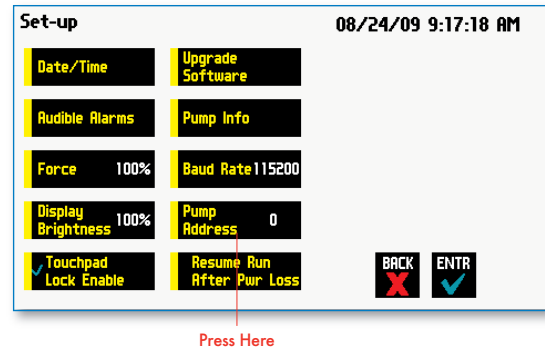
Up to 100 pumps can be linked together. Each pump can be given a unique address.

In the Set-up Screen press **Pump Address**.
Valid entries are 00 to 99.

Enter the pump address.

Accept the changes made by pressing **ENTR** and go back to the Set-up screen.

Cancel the changes made by pressing **BACK** and go back to the Set-up screen.



Setting the RS 232 Baud Rate

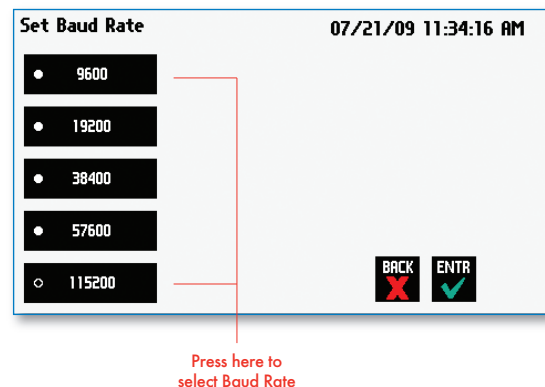
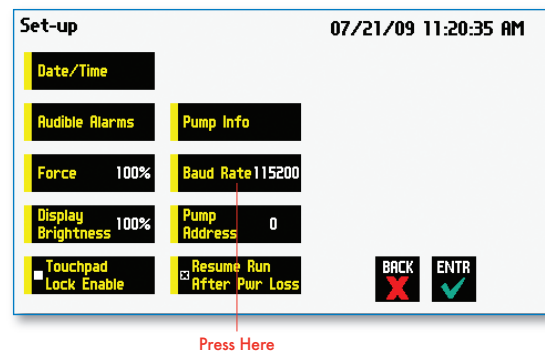
The Cole-Parmer® Touch Screen Series has an RS 232 port for serial communications.

To adjust the baud rate from 9600 to 115,200: in the Set-up screen press **Baud Rate**.

Select the Baud Rate.

Accept the changes made by pressing **ENTR** and go back to the Set-up screen.

Cancel the changes made by pressing **BACK** and go back to the Set-up screen.



Resume Run on Power Loss

If Resume Run on Power Loss is selected the pump will power up running when power resumes.



NOTE: This feature is only valid in Infuse Only or Withdraw Only modes. This feature will only execute if the pump is running at the time when power is lost, and was running with no target set.

In the Set-up screen, select **Resume Run After Power Loss**



Accept the changes made and go back to the previous screen before Set-up mode was selected.



Cancel and go back to the previous screen before Set-up mode was selected.

Set-up 08/24/09 9:17:18 AM

Date/Time	Upgrade Software
Audible Alarms	Pump Info
Force 100%	Baud Rate 115200
Display Brightness 100%	Pump Address 0
Touchpad Lock Enable	Resume Run After Pwr Loss

BACK ENTR

Press Here



NOTE: Elapsed Timers and Total Infused Volume will be reset when power resumes. Counters and timers information prior to loss of power are lost

Setting and Viewing Pump Information

The information on the pumps parameters are displayed in the Pump Information Screen.

In the Set-up screen, select **Pump Information**

The current information on the pump will be displayed.

- The Model type
- The Software Version
- The # of program steps used out of the # of program steps available
- The Calibration Due Date
- The Maintenance Due
- The Current Pump Address

Press **Calibration Reset** to reset the calibration date (one year from current date).

Confirm the Calibration Reset by pressing the **Calibration Reset** again.

Press **Maintenance Reset** to reset the # of hours for maintenance. This will reset it to 100 hours.

Confirm the Maintenance Reset by pressing the **Maintenance Reset** again.

Set-up 08/24/09 9:17:18 AM

Date/Time	Upgrade Software
Audible Alarms	Pump Info
Force 100%	Baud Rate 115200
Display Brightness 100%	Pump Address 0
Touchpad Lock Enable	Resume Run After Pwr Loss

BACK ENTR

Press Here

Pump Information 10/30/09 02:16:39 AM

Model Type	Programmable
Software Version	0.0.7e
Program Steps Available	756/807
Calibration Due Date	10/01/10
Maintenance Due	85 hours
Pump Address	0

Calibration Reset Maintenance Reset

BACK ENTR

Press here to reset the calibration date

Press here to reset the maintenance date



NOTE: The calibration date is always shown in month/days/year.

Accept the changes made by pressing and go back to the Set-up screen.



Cancel the changes made by pressing and go back to the Set-up screen.

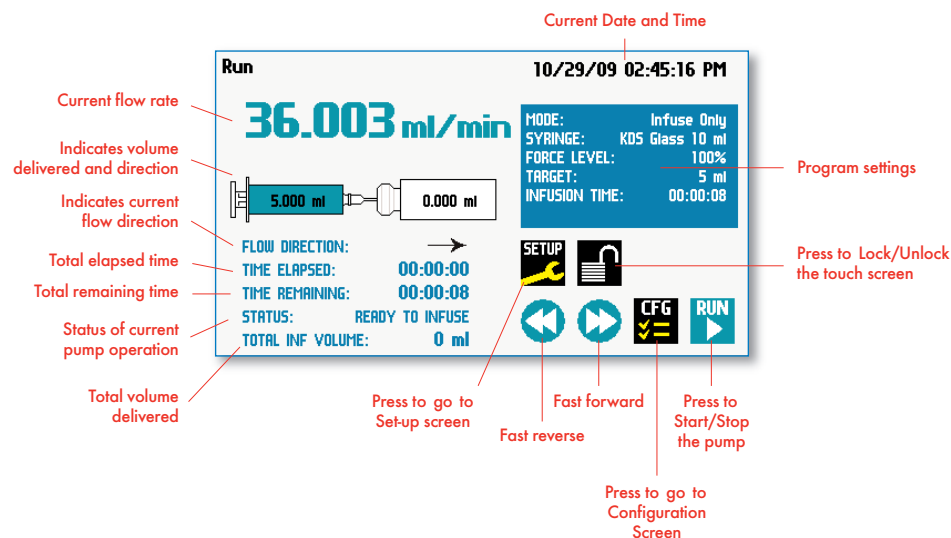


NOTE: Free RTOS v5.4.2 (www.freertos.com)

Operating the Pump

The Cole-Parmer® Touch Screen Series is easy to operate because all the pump parameters are shown on the RUN Screen. Depending on the mode of the pump, different parameters will be shown on the display.

The Cole-Parmer® Touch Screen Series Run Screen is shown below.



Running the Pump

Simply press the **RUN** button to start the pump.



When the pump is running the **STOP** button is shown.

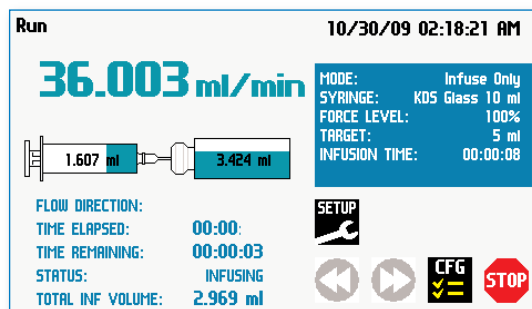
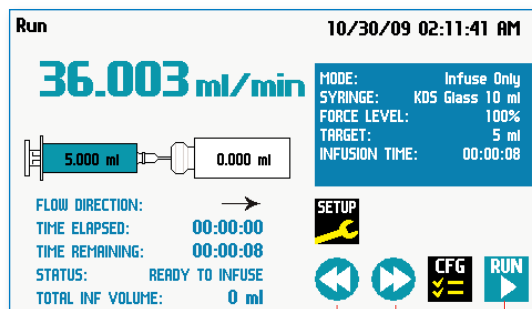


Press the **STOP** button to pause the program.




The fast forward and fast reverse key are not active when the pump is running.

The Cole-Parmer® Touch Screen Series can also be run in a vertical position. All screens will rotate automatically when the pump is placed on its side.



Custom Program Workflow

The following flow diagram describes the use of an Custom Program:

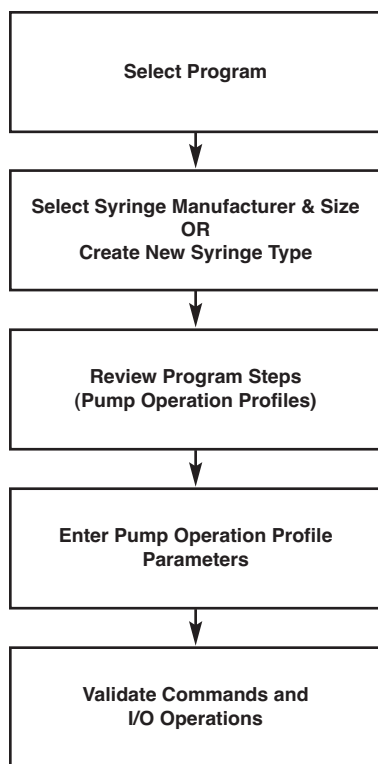

The software provides the user with a variety of prompts and warnings to help guide the setup of programs. For example, if a pump Program incorporates steps that exceed the capacity of the defined syringe, the system will prompt you that parameters are out-of-range. The system also prompts users if required values are not entered.



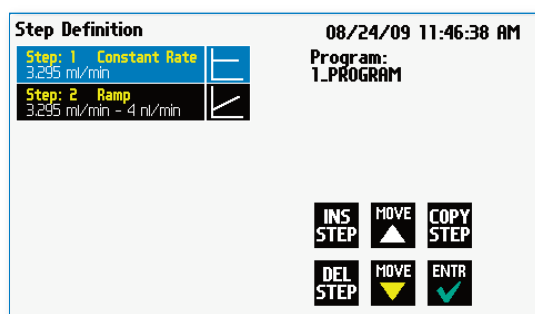
NOTE: The pump will check the validity of parameters entered for each Step in the Program, but will not review the Program in its entirety to verify that it is logical and can be fully executed

STEP-BY-STEP OPERATION EXAMPLE: USING A CUSTOM PROGRAM

The following steps describe operation using a Custom Program.

In the following description, we are assuming the use of a Program that is already saved on the Cole-Parmer® Touch Screen Series. Note that if you make any changes to the Program's operating parameters and choose Accept on the page where the changes are made, the changes automatically become part of the stored Program for future use.

1. From the Program Configuration screen, the Program you wish to use should be displayed in the Program Select button. If the Program is not displayed, select the correct Program using the procedure outlined in Running The Pump.
2. Choose the Syringe Select button to enter the Syringe Selection screen.
3. Follow the instructions in the Syringe Selection section of this Guide to complete the selection of a syringe from the built-in syringe definition list, or to create a new syringe definition. When finished, choose the Accept button to save your changes, or choose the Cancel button to cancel your changes, and return to the Program Configuration screen.
4. Review the Pump Control Profile Steps described in the Step Definition button. Note that the software will display the first three Steps in the Program in the main button area. If the Program uses more than three Steps, the last Step in the Program will be displayed in the box that is still part of the button beneath the Step Definition button.
5. To review the Program Steps in detail or to update/modify any of the stored Pump Control Profile Parameters, choose the Step Definition button to enter the Program Steps Summary Screen. The name of the Program is listed at the right of the screen, and the complete list of Steps is listed in a series of buttons is on the left of the screen.
6. To review details for any individual Step in the Program, select the Step you wish to review, then press the Step button again to enter the Step Definition screen for that Step.
7. To edit the parameters for any Pump Control Profile, follow the instructions supplied in the applicable Pump Operation Profile description contained in the Setting Program Parameters section later in this Guide.
8. To return to the Program Steps Summary screen when finished reviewing or editing the Step parameter details, choose the Accept button to save your changes or choose the Cancel button to cancel any changes you made.
9. When you are finished reviewing the operating details for all Steps, choose the Back button to return to the Program Configuration screen.
10. To preview the Program, choose the Preview Program button, or to run the Program, choose the Program Run button.



Step Definition Screen showing the list of Steps currently defined for the Program

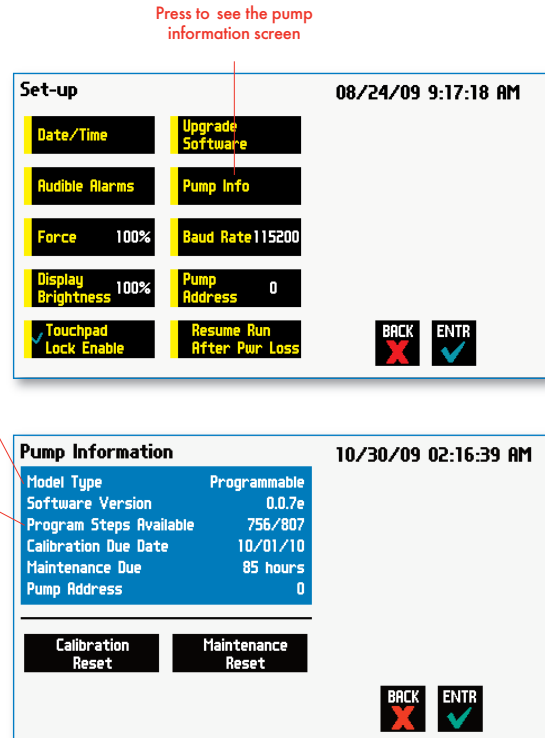
Configuring Custom Programs

The Cole-Parmer® Touch Screen Series Programmable Models offer the flexibility of configuring custom programs including multistep programs. The programs and steps can also be triggered to start and stop with events and external I/O, and user defined inputs. The maximum number of steps is 800 for all the programs.

To determine if your pump is programmable go to the pump information screen. The Model Type should say Programmable.

Press 


The available number of program steps is shown in the pump information screen.



Configuring a New Custom Program


To display the **Mode Selection** screen press Mode.

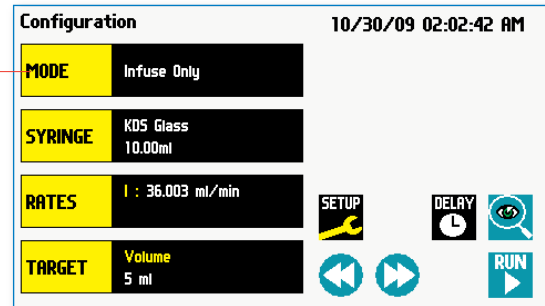
Select Mode





The Mode Selection screen will be displayed

Select **New Program**

View programs up a page 

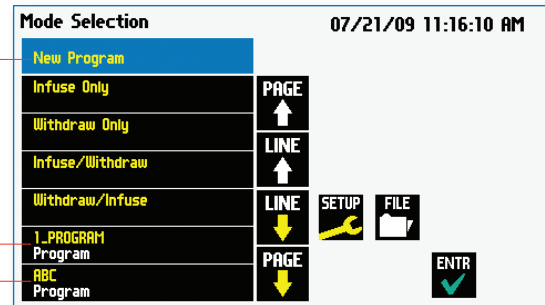
View the programs down a page 

Views up the program 

Views down the program 

Press

Custom Program



Configuring Custom Programs

The first thing to do on a custom program is to name the program.

Enter any 15 character alphanumeric name

Accept the new name



Go back to the [Program Selection](#) screen



Select the Syringe manufacturer and Size
(see page 16)

To define the steps for the custom program
press [Step Definition](#)

Below are some predefined steps including:

Constant Rate- a single continuous rate of flow

Ramped Rate - a linearly changing flow rate between
a starting and ending flow rate

Stepped Rate- an evenly stepped changing flow
rate between a starting and ending flow rate

Pulsed Rate- an oscillating flow rate varying between
a starting and a target rate

In addition to these you can repeat, go back or jump
forward in a program as well as send signals to external
devices or react to signals received from external devices.
This package of pumping profiles and commands
provides users with an extremely broad array of
applications that can be easily configured with complex
external programming.

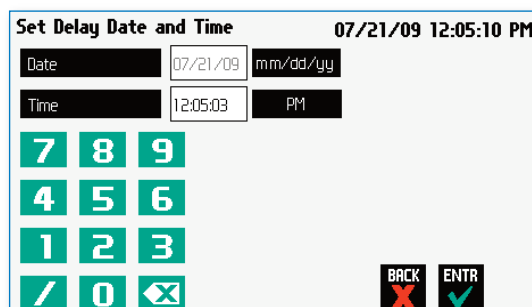
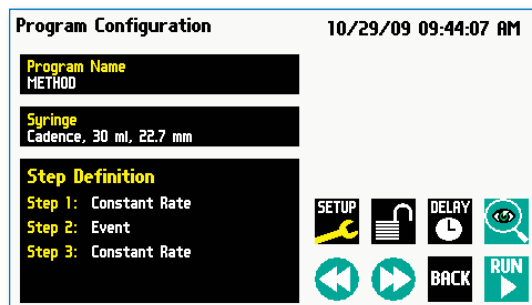
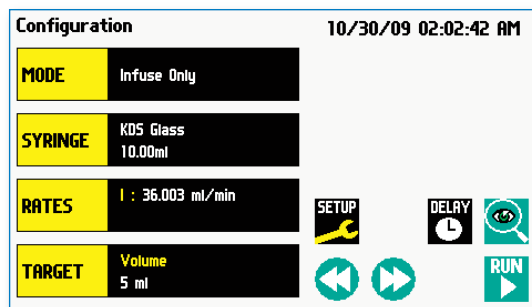
Programs can be edited and saved under a new name,
simply edited to fit another application. Steps can be
copied, moved, inserted or deleted.

Setting A Time Delay


You can set the Cole-Parmer® Touch Screen Series to run at a specific date and time in the future. This delayed run function allows you to set up all the parameters for operation, then have the pump start unattended at the required time.

SETTING A TIME DELAY FOR PROGRAM OPERATION:

1. From the Program Configuration screen, choose the Set Delay Parameters button to enter the Set Delay Date and Time screen.
2. To enter the date, select the Date button and type the desired date using the onscreen number pad. Select the Date Format button to toggle between Month/Day/Year (m/d/y), Day/Month/Year (d/m/y) and Year/Month/Day (y/m/d) formats. Enter your values in two-digit format for each value. For example, April 3, 2009 would be entered as 04/03/09 in m/d/y format.
3. To enter the time, select the Time button and type the desired time using the onscreen number pad. Select the Time Period button to toggle between AM/PM and 24-hours settings. Time should be entered using hours:minutes:seconds (hh:mm:ss).
4. When the desired values are entered, choose the Accept button to save your Changes and return to the Program Configuration screen. When a Time Delay value is active, the Time Delay Parameters button is displayed in Red on the screen to alert the user that a delay has been set.
5. Press the Run button to activate the Delay.




Set Delay Date and Time screen

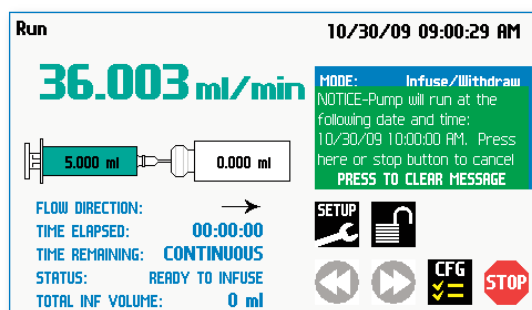

Quick Start run screen with a Time Delay Set

Creating a Custom New Program

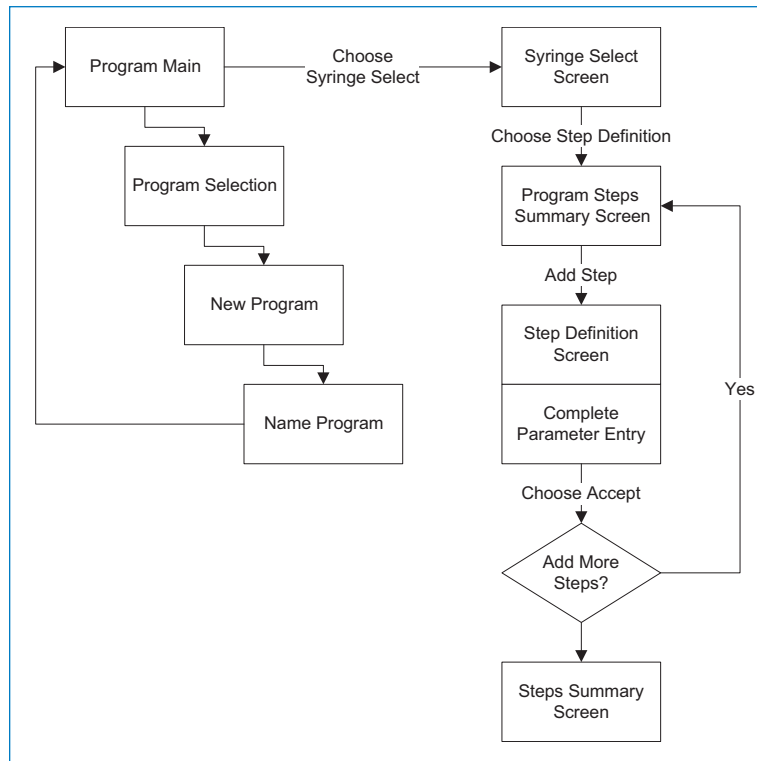If none of the Programs stored on the Cole-Parmer® Touch Screen Series match your requirements, you can create a new Program and store the Program for use. The following flow diagram outlines the basic steps you will follow to create a new Program:



New User-Defined Program creation workflow

TO CREATE A NEW USER-DEFINED PROGRAM

1. From the Program Configuration or Programs screen, choose the Program Select button to enter the Program Selection screen.
2. Select New Program then choose Accept or press New Program again.
3. On the Program Name screen, type the name of your Program using the onscreen keyboard. Choose the Symbols/Numbers button **.?123** to display the list of available numbers and symbols that can be included in the Program Name. To switch the keyboard back to letters mode, choose the Letters button (ABC).




New Program Name Entry screen showing letter keys.
To view numbers/symbols, choose the **.?123** use button in the lower left corner

TO CREATE A NEW USER-DEFINED PROGRAM (CONTINUED)

4. When you have entered the desired name for your Program, choose the Accept button to go to the Program Configuration screen. Your new Program will now be listed in the Program Select button.
5. Choose the Syringe Select button and complete the description of the syringe by following the procedure outlined previously in the Syringe Selection section of this Guide.
6. To begin creating Pump Control Profile Steps, choose the Step Definition button to enter the Program Steps Summary screen.
7. To add a Step to your Program, choose the Add Step button to access the Program Steps Menu screen. The available options include:

- | | |
|------------------|---------------------|
| a. Constant Rate | d. Pulse |
| b. Ramp | e. Advanced Options |
| c. Stepped | |

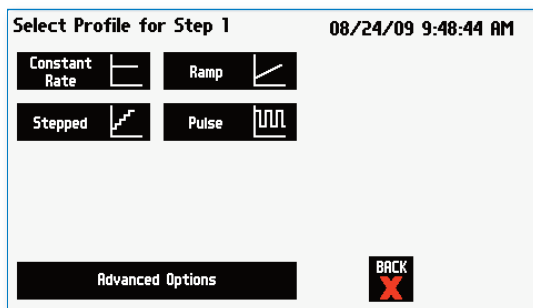

Program Steps Profile Menu screen

8. Choose the Program Step type that you want to add to the Program. If you choose a pumping profile (e.g. Constant Rate, Ramp, etc.) you will be taken to the Step Definition screen for the selected Program. If you choose Advanced Options, you will be taken to the Select Advanced Options screen.
9. Complete the parameter entry applicable to the selected Pump Operation Profile, follow the detailed instructions supplied in the Pump Operation Profiles (Existing and User-Defined Programs) section later in this Guide. If you selected Advanced Options, follow the instructions supplied in the Incorporating Advanced Options Into a Program section later in this manual.



NOTE: If you do not complete the required parameter entry for a Step, subsequent Step additions may be ignored. You should enter all required parameters for a Step when creating it to ensure that the application software can assemble your Program without errors.

10. When you have completed all parameter entries on the Step Definition screen or completed the required entries for Advanced Options, choose the 'BACK' button to save your changes and return to the Program Steps Summary screen.



11. Repeat items 7-10 in this procedure to include additional Steps in the Program.



NOTE: When adding a Step to a Program, the Step is added after the currently selected Step (i.e. the Step that is highlighted in blue)

12. When you are done adding all required Steps to the Program, choose the 'BACK' button to return to the Program Main screen.

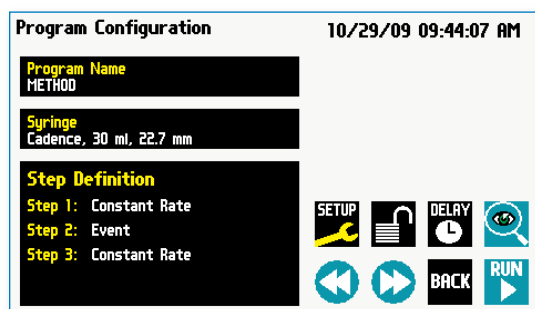


Modifying an Existing User-Defined Program

Once you have created a Program or have reviewed an Existing Program, you may have a need to modify the structure of the Program (i.e. the types of Steps and Commands included, not just changes in the operating parameters for a specific Step. The following procedures will help you to efficiently modify a stored Program.

CHANGING THE ORDER OF STEPS IN A PROGRAM


1. From the Program Configuration screen, choose the Step Definition button to enter the Program Steps Summary screen.



2. Highlight (select with a single touch) the Step for which you want to change the order.
3. Choose the Move Step Up or Move Step Down button, as required, to move the Step into the position that you want it executed.




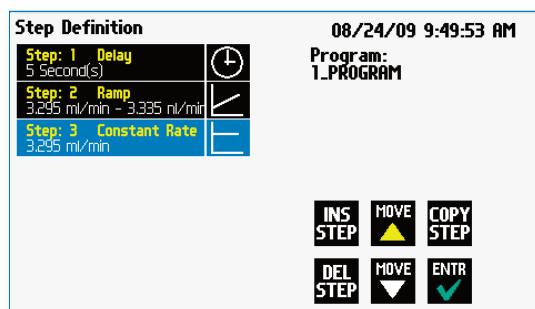NOTE: If there are Program Commands (e.g. Go To) that link to the Step, moving Steps can result in unexpected results as the Command may be linked to the Step number. You should always recheck all Program Commands after adding, deleting or moving Steps in a previously created Program.



DELETING A STEP FROM A PROGRAM

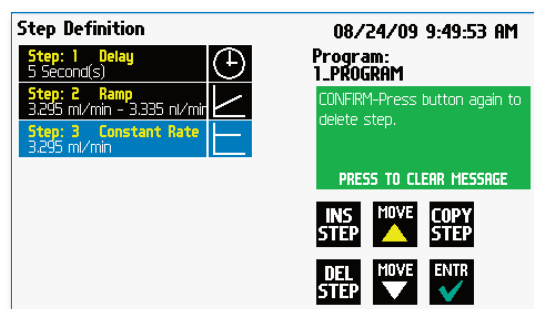1. From the Program Configuration screen, choose the Step Definition button to enter the Program Steps Summary screen.
2. Highlight (select with a single touch) the Step that you want to delete.
3. Choose the Delete Step button. The application will display a warning message in the Message Area asking you to confirm your deletion request.



Delete Program Step showing confirmation message on the screen

4. To confirm the deletion, choose the Delete Step button again. The Step is removed and you are returned to the Program Steps Summary screen.

Defining Program Step Parameters

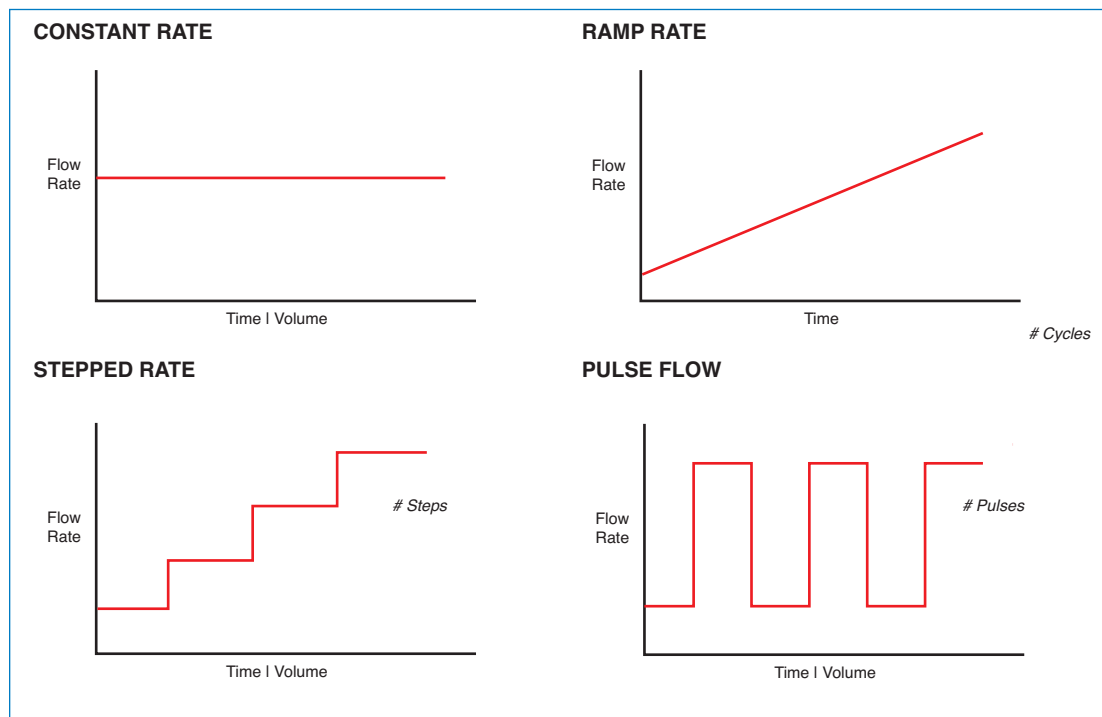
This section of the User Guide describes each of the available functions available for use in Cole-Parmer® Touch Screen Series Existing and User-Defined Programs. The information in this section is designed to help you set the individual parameter values in each of the Program Step options and to properly use the Advanced Options commands.

The information that follows is organized into two sections:

- Pump Operation Profiles
- Advanced Options (Commands and I/O)

PUMP OPERATION PROFILES


The Cole-Parmer® Touch Screen Series Existing Programs and all User-Defined Programs use Pump Operation Profiles as part of the Programs controlling the pump. These profiles are combined with a syringe selection and any Commands or Input/Output triggers to create a complete Program. The following sections describe each profile along with instructions for setting up the profile. The following Pump Operation Profiles are available on the Cole-Parmer® Touch Screen Series:



Pump Operation Profiles available on the Cole-Parmer® Touch Screen Series

Constant Rate Profile

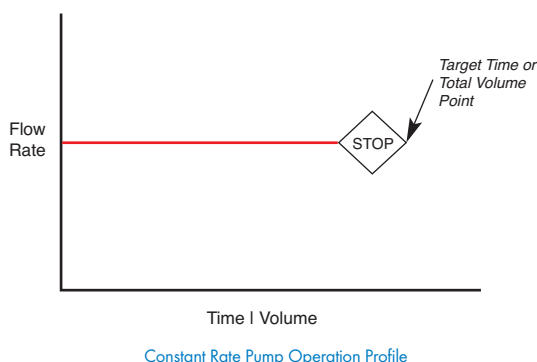The Constant Rate profile is used to run the pump at a defined flow rate and for a specific pumping mode (infuse or withdraw) up to a target volume or for a defined time period.



When creating or editing a Constant Rate profile, you will define three primary operating parameters:

- The Flow Mode defining the direction of flow: infuse or withdraw
- The Flow Rate to define the rate of uptake or dispensing
- The Target Volume or Run Time that determine how long the pump will operate at the set rate

When creating or editing a Constant Rate profile, you can set the rate in ml, μ l, nl or pl units of volume using sec, min or hr as the units of time. The system automatically calculates a minimum and maximum flow rate based on the syringe selected or custom-defined. The min/max values are made available as button selection options on the Set Infuse and Withdraw Rate screens.

The screenshot shows the 'Step # 4, Constant Rate' setup screen. It includes a timestamp '07/21/09 12:13:28 PM'. The 'Direction' is set to 'Infuse' with a right-pointing arrow. The 'Rate' is set to '0 ml/min' and the 'Volume' is set to '0 ml'. At the bottom, there are 'BACK' and 'ENTR' buttons.

Constant Rate Profile Setup Screen

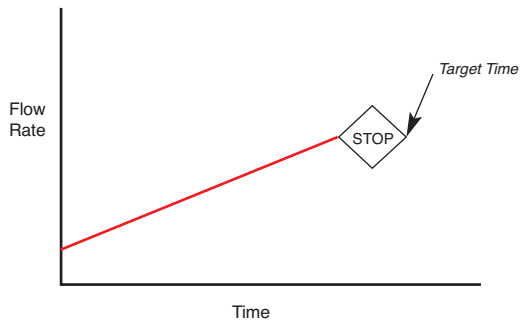
TO SET UP/EDIT A CONSTANT RATE PROFILE:

1. From the Step Definition screen for a Constant Rate profile, choose the Direction button to toggle the Flow Mode between Infuse and Withdraw operation.
2. Choose the **Rate** button to enter the Set Rate screen.
3. Enter the flow rate value desired using the keypad or select the **Minimum Rate** or **Maximum Rate** button.
4. If the Units of Rate are not correct for your application, choose the **Units** button.
5. On the Units Selection screen, select the unit of measure (ml, μ l, nl or pl) and the unit of time (sec, min, hr) as appropriate, then choose 'ENTR' to confirm your selections.
6. When all parameters are correctly set on the Set Rate screen, choose 'ENTR' to confirm your settings and return to Constant Rate Definition screen.
7. Choose the **Volume** button to enter the Set Volume or Time screen.
8. On the Set Volume or Time screen, select the **Volume or Time** button, as desired.
9. Enter the desired value using the keypad, or select the **Maximum Volume** button.
10. If the Units of Volume/Time are not correct for your application, choose the **Units** button to the right of the parameter entry box. For the Units of Volume, this choice takes you to the Units Selection screen. For the Units of Time option, this choice will toggle the button between sec (seconds) and h:m:s (hours:minutes:seconds). To use the h:m:s option, enter your time value in hours, minutes and seconds, separating each value with a colon.
11. If you are changing the Units of Volume, on the Units Selection screen, select the unit of measure (ml, μ l, nl or pl) as appropriate, then choose 'ENTR' to confirm your selection.
12. When all parameters are correctly set on the Set Target Volume or Time screen, choose 'ENTR' to confirm your settings and return to the Constant Rate Definition screen.
13. Choose 'ENTR' again to return to the Step Definition screen.



Ramp Profile

The Ramp Rate profile is used to create a linearly increasing or decreasing flow rate for a single pumping mode (infuse or withdraw) over the course of a defined time period.



Ramp Rate Pump Operation Profile

When creating or editing a Ramp Rate profile, you will define three primary operating parameters:

- The Flow Mode defining the direction of flow: *infuse or withdraw*
- The Flow Rate at the beginning and end of the pumping cycle
- The Target Time that determines the total length of time over which the ramp will be calculated

During a Ramp Rate profile, the Cole-Parmer® Touch Screen Series varies the flow rate linearly between the starting rate and ending rate over the full time period entered. A Ramp Rate profile can be applied to an increasing or decreasing rate application.

When creating or editing a Ramp Rate profile, you can set the rate in ml, µl, nl or pl units of volume using sec or h:m:s (hours:minutes:seconds) as the units of time. The system automatically calculates a minimum and maximum flow rate based on the syringe selected or custom-defined. The min/max values are made available as button selection options on the Set Infuse and Withdraw Rate screens.

TO SET UP/EDIT A RAMP RATE PROFILE:

1. From the Step Definition screen for a Ramp Rate profile, choose the Direction button to toggle the Flow Mode between Infuse and Withdraw operation.

Step # 4, Ramp 07/21/09 12:16:39 PM

Direction Infuse →

Rates
Start: 0 ml/min End: 0 ml/min

Time
0 second(s)

BACK ENTR

Step Definition screen for a Ramp Rate profile

2. Choose the Set **Start & End Rate** button to enter the Set Start & End Rates screen.

Set Start and End Rate 08/24/09 10:07:09 AM

Start Rate 0 ml/min

End Rate 0 ml/min

7 8 9 ↑ MAX 3.452 ml/min

4 5 6

1 2 3 3.325 nl/min

0 0 MIN ↓

BACK ENTR

Set Start and End Rates screen

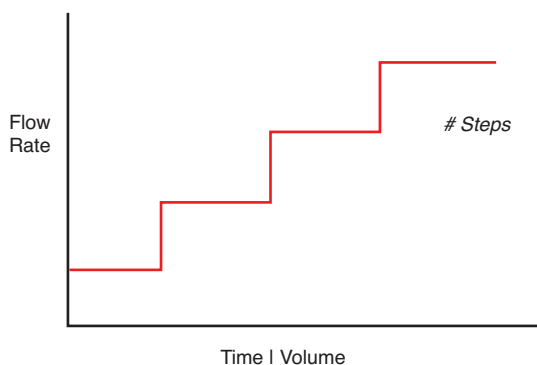
3. Select the **Start Rate** button (it should be selected by default as indicated by the active data entry box to the right of the button). Enter the flow rate value desired using the keypad or select the **Minimum Rate** or **Maximum Rate** button.
4. If the Units of Rate are not correct for your application, choose the **Units** button.
5. On the Units Selection screen, select the unit of measure (ml, µl, nl or pl) and the unit of time (sec, min, hr) as appropriate, then choose **ENTR** to confirm your selections.
6. Repeat Steps 3 to 5 for the End Rate by selecting the **End Rate** button and continuing data entry.
7. When all parameters are correctly set on the Set Rate screen, choose **ENTR** to confirm your settings and return to Ramp Definition screen.
8. Choose the Time button to enter the Set Time screen.
9. Enter the desired Time value using the keypad.
10. If the Units of Time are not correct for your application, select the **Units** button to toggle between sec (seconds) and h:m:s (hours:minutes:seconds). To use the h:m:s option, enter your time value in hours, minutes and seconds, separating each value with a colon.
11. When all parameters are correctly set on the Set Target Time screen, choose Accept to confirm your settings and return to the Ramp Definition screen.
12. Choose Accept again to return to the Step Definition screen.



NOTE: Ramps ≤ 2 seconds from min to max rates may cause motor stalls.

Stepped Profile

The Stepped Rate profile is used to create an increasing or decreasing flow rate that varies in a series of equal steps between a starting and ending flow rate.



Stepped Rate Pump Operation Profile

When creating or editing a Stepped Rate profile, you will define four primary operating parameters:

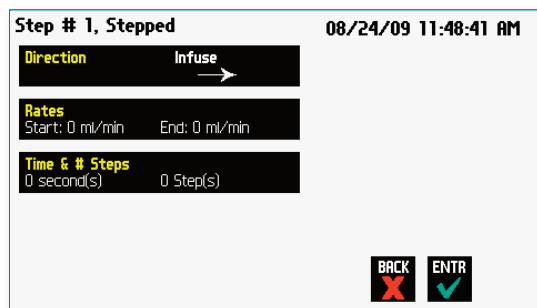
- The Flow Mode defining the direction of flow:
infuse or withdraw
- The Flow Rate at the beginning and end of the pumping cycle
- The Target Time that determines the total length of time over which the stepped profile will be calculated
- The Number of Steps that will be created in the profile

During a Stepped Rate profile, the Cole-Parmer® Touch Screen Series calculates a series of equal rate increases that will vary the flow from the starting value to the ending value. You enter the Number of Steps you want to use. The Target Time defines the total length of pumping time, with the Target Time value divided equally into each stepped rate. A Stepped Rate profile can be applied to an increasing or decreasing rate application.

When creating or editing a Stepped Rate profile, you can set the rate in ml, µl, nl or pl units of volume using sec or h:m:s (hours:minutes:seconds) as the units of time. The system automatically calculates a minimum and maximum flow rate based on the syringe selected or custom-defined. The min/max values are made available as button selection options on the Set screen.


TO SET UP/EDIT A STEPPED RATE PROFILE:

- From the Stepped Definition screen for a Stepped Rate profile, choose the **Direction** button to toggle the Flow Mode between Infuse and Withdraw operation.



Step Definition screen for a Stepped Profile

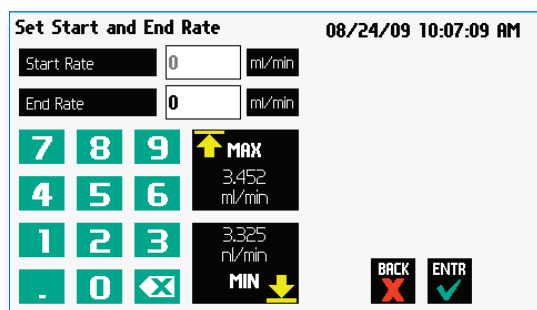- Choose the **Rates** button to enter the Set Start and End Rates screen.




Set Start and End Rate screen

- Select the **Start Rate** button (it should be selected by default as indicated by the active data entry box to the right of the button). Enter the flow rate value desired using the keypad or select the **Minimum Rate** or **Maximum Rate** button.
- If the Units of Rate are not correct for your application, choose the **Units** button.
- On the Units Selection screen, select the unit of measure (ml, µl, nl or pl) and the unit of time (sec, min, hr) as appropriate, then choose Accept to confirm your selections.
- Repeat Steps 3 to 5 for the End Rate by selecting the **End Rate** button and continuing data entry.
- When all parameters are correctly set on the Set Start and End Rate screen, choose Accept to confirm your settings and return to Stepped Definition screen.
- Choose the **Time** and **# Steps** button to enter the Set Time & Steps screen.

- Select the **Total Time** button (it should be selected by default as indicated by the active data entry box to the right of the button). Enter the desired value using the keypad.
- If the Units of Time are not correct for your application, select the **Units** button to toggle between sec (seconds) and h:m:s (hours:minutes:seconds). To use the h:m:s option, enter your time value in hours, minutes and seconds, separating each value with a colon.
- Select the **Number of Steps** button and enter the desired value using the keypad.
- When all parameters are correctly set on the Set Time and Steps screen, choose Accept to confirm your settings and return to the Stepped Definition screen.

Pulse Profile

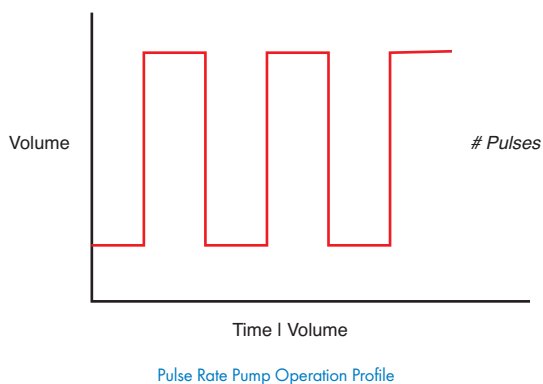The Pulse Rate profile is used to create an oscillating infusion or withdrawal profile that varies between a starting flow rate and a second rate that can be either or higher or lower than the start rate. The flow rate is constant within each segment of the pulsed flow profile and the user defines the total number of repetitions (pulses) for the oscillating pattern. Each segment of the profile must be set to run for either a specific length of time or to deliver a set volume.



When creating or editing a Pulse Rate profile, you will define four primary operating parameters:

- The Flow Mode defining the direction of flow:
infuse or withdraw
- The Flow Rate to be used at each level of the pulsed profile within the pumping cycle
- The length of Time that each segment of the profile will run or the specific Volume to be infused/withdrawn during each segment of the profile
- The Number of Pulses (repetitions) that will be executed during the pumping cycle

During a Pulse Rate profile, the Cole-Parmer® Touch Screen Series calculates a series of equal length pumping stages at two user-defined flow rates: a starting rate and a “pulse” rate. You enter the Number of Pulses (repetitions) that you want to use, along with how you want the system to calculate the length of each stage – by time or by volume pumped. A Pulse Rate profile can be applied to pulses that are higher in rate or lower in rate than the starting flow rate.

When creating or editing a Pulse Rate profile, you can set the rate in ml, µl, nl or pl units of volume using sec, min or hr as the units of time. The system automatically calculates a minimum and maximum flow rate based on the syringe selected or custom-defined. The min/max values are made available as button selection options on the Set Pulse Rates screens.

TO SET UP/EDIT A PULSE RATE PROFILE:

- From the Pulse **Definition** screen for a Pulse Rate profile, choose the **Direction** button to toggle the Flow Mode between Infuse and Withdraw operation.

Step # 4, Pulse 07/21/09 12:28:51 PM

Direction Infuse →

Rates
Rate 1: 0 ml/min Rate 2: 0 ml/min

Volume
0 ml 0 ml

○ Volume ● Time

of Pulses
0 Pulse(s)

BACK ENTR

Step Definition screen for a Pulse Profile
with Volume mode selected

Set Pulse Times 08/24/09 10:09:02 AM

Time 1 0 SECS

Time 2 0 SECS

7 8 9
4 5 6
1 2 3
- 0

BACK ENTR

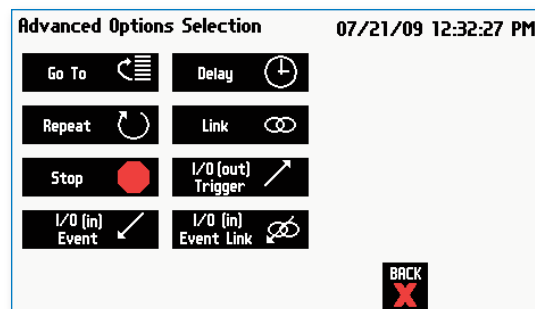
Set Pulse Times screen used to define the length of time
the pump will operate at each pulse step level

- Choose the **Rates** button to enter the Set Pulse Rates screen.
- Select the **Rate 1** button (it should be selected by default as indicated by the active data entry box to the right of the button). Enter the flow rate value desired using the keypad or select the **Minimum Rate** or **Maximum Rate** button.
- If the Units of Rate are not correct for your application, choose the **Units** button.
- On the Units Selection screen, select the unit of measure (ml, µl, nl or pl) and the unit of time (sec, min, hr) as appropriate, then choose Accept to confirm your selections.
- Repeat Steps 3 to 5 for the second flow rate by selecting the **Rate 2** button and continuing data entry.
- When all parameters are correctly set on the Set Pulse Rates screen, choose Accept to confirm your settings and return to Pulse Definition screen.
- Select the **Volume** or **Time** button with indicator to determine how the pump cycle length will be managed. Volume calculates each cycle stage by pumping for the length of time needed to pump the specified volume. Time manages each cycle stage by pumping for the amount of time specified.
- Depending on whether you have selected the **Volume** or **Time** button in the preceding step, choose the **Time** or the **Volume** button to enter the Set Pulse Times or the Set Pulse Volumes screens.
- If you selected Set Time, continue with Step 11. If you selected Set Volume, proceed to Step 14.
- On the Set Pulse Times screen, select the **Time 1** button (it should be selected by default as indicated by the active data entry box to the right of the button). Enter the desired value using the keypad. Repeat for Time 2.
- If the Units of Time associated with Time 1 or Time 2 are not correct for your application, select the **Units** button associated with either parameter to toggle between sec (seconds) and h:m:s (hours:minutes:seconds). To use the h:m:s option, enter your time value in hours, minutes and seconds, separating each value with a colon.
- When both Time parameters are set to the correct values, choose the Accept button to return to the Pulse Definition screen. Continue with Step 17.
- On the Set Pulse Volumes screen, select the **Volume 1** button (it should be selected by default as indicated by the active data entry box to the right of the button). Enter the desired value using the keypad. Repeat for Volume 2.
- If the Units of Volume associated with **Volume 1** or **Volume 2** are not correct for your application, choose the **Units** button associated with either parameter to enter the Units selection screen. Select the unit of measure button (ml, µl, nl, pl) applicable to your application, then choose the Accept button to return to the Set Pulse Volumes screen.
- When both Volume parameters are set to the correct values, choose the Accept button to return to the Pulse Definition screen.
- Select the **# of Pulses** button to enter the Set # Pulses screen.
- On the Set # Pulses screen, enter the desired value using the keypad, then choose the Accept button.

Incorporating Advanced Options (Commands and I/O) Into A Program

Cole-Parmer® Touch Screen Series Programs can incorporate a variety of commands that expand the versatility of the system by permitting convenient looping, chaining of events, sending signals to external devices, or accepting inbound signals from external devices. These Advanced Options controls are described in this section of the Guide.

To access the Advanced Options, choose Advanced Options from the Select Profile screen.



Program Step Profile Selection Advanced Options screen showing available command options

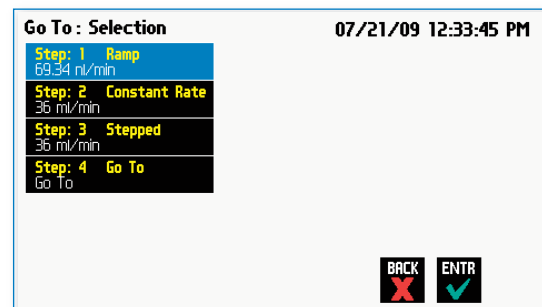
Go To Step

The Go To Command is used to instruct the Cole-Parmer® Touch Screen Series to jump to a particular Step in a Program. This can be used to help you loop back or jump forward in a Program protocol and can be used in conjunction with I/O Events to make the pump respond to a trigger from an external device.

When defining a Go To jump in a Program, the Command targets the Step sequence position and not the function that currently occupies that position. For example, if you create a Go To jump to Step 4 in a Program (currently a Ramp Pump Operation Profile) and subsequently move a Constant Rate profile into the Step 4 position, the Go To jump will then connect to the Constant Rate operation and not to the Ramp operation.

TO INSERT A GO TO JUMP:

1. From the Program Main screen, choose Step Definition to enter the Step Definition screen.
2. Highlight (select once) the Program Step that you want the Go To jump to follow.
3. Choose the Insert Step button to access the Program Steps Menu screen.
4. Choose Advanced Options to enter the Set Advanced Options screen.
5. Choose Go To to insert the Go To Step and view an updated list of all current Steps in the Program.



Advanced Options Go To Step screen. Choosing the Go To Command brings you to a listing of the available steps in the Program

6. Double tap the Step that you want to jump to with the Go To Command. The system automatically returns to the Step Definition screen which now shows the Go To Step inserted in the Program.

Delay

The Delay Command is used to insert a time delay in a procedure. You can program a time delay between 0.2 seconds up to 99:99:99 (99 hours, 99 minutes, 99 seconds).

TO INSERT A TIME DELAY IN A PROGRAM:

1. From the Program Main screen, choose Step Definition to enter the Step Definition screen.
2. Highlight (select once) the Program Step that you want the time delay to follow.
3. Choose the Insert Step button to access the Program Steps Menu screen.
4. Choose Advanced Options to enter the Set Advanced Options screen.
5. Choose Delay to enter the Set Delay Time screen.

Advanced Options Set Delay Time screen

6. Enter the desired Time using the keypad.
7. If the Units of Time associated with Time are not correct for your application, select the Units button to toggle between sec (seconds) [0.2 seconds minimum allowable value] and h:m:s (hours:minutes:seconds). To use the h:m:s option, enter your time value in hours, minutes and seconds, separating each value with a colon. You can enter a value up to 99:99:99.
8. When your Time entry is correct, choose the Accept button to insert the Delay Step in the Program and return to the Step Definition screen.

Repeat Step

The Repeat Command is used to repeat Steps in a procedure a defined number of times. To use this Command, you will tell the application what Steps you want to Repeat and how many times you want them repeated. The Repeat Step command creates a "jump" back to a specific Step in a Program that executes the defined Step plus all following Steps up to the Repeat command (i.e. it creates a Program loop).

When defining a Repeat jump in a Program, the Command targets the Step sequence position and not the function that currently occupies that position. For example, if you create a Repeat jump to Step 4 in a Program (currently a Ramp Pump Operation Profile) and subsequently move a Constant Rate profile into the Step 4 position, the Repeat jump will then connect to the Constant Rate operation and not to the Ramp operation.

TO INSERT A STEP REPEAT SEQUENCE IN A PROGRAM:

1. From the Program Main screen, choose **Step Definition** to enter the Step Definition screen.
2. Highlight (select once) the Program Step that you want the repeat step to follow.
3. Choose the **'INS STEP'** button to access the Program Steps Menu screen.
4. Choose Advanced Options to enter the Set Advanced Options screen.
5. Choose Repeat to insert the Repeat Step function into the Program.
6. Select the Step(s) you would like to repeat from by selecting the Step button twice or selecting once and then choosing Accept.

INS STEP

Advanced Options Repeat Step screen. Choosing Repeat brings you to a listing of the available steps in the Program

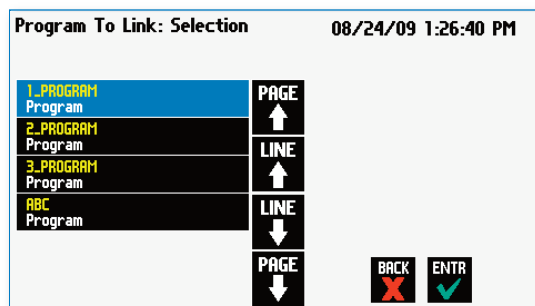
7. Enter the number of repeats using the keypad then choose the Accept button to save your change and return to the Step Definition screen.



NOTE: This command creates a loopback to the designated Step in a Program, executing the selected Step, plus all subsequent Steps up to the point of insertion of the Repeat command. For example, if the Repeat command is set to repeat Step 2 one time, the system will loop back to Step 2 and execute both Steps 2 and 3 one additional time.

Link Program

The Link Command is used to insert a link to another stored Program into the operating sequence of a Program. You can Link a Program any place in the Program sequence.



Advanced Options Link Command selection screen

TO LINK TO A STORED PROGRAM:

1. From the Program Main screen, choose Step Definition to enter the Step Definition screen.
2. Highlight (select once) the Program Step that you want the Link Command to follow.
3. Choose the 'INS STEP' button to access the Program Steps Menu screen.
4. Choose **Advanced Options** to enter the Set Advanced Options screen.
5. Choose Link to enter the Link Selection screen and display a list of the saved Programs available.
6. Double tap the Program that you want to link to in the current Program. If the desired Program is not listed on the screen, use the 'LINE UP' 'LINE DOWN' or 'PAGE UP' 'PAGE DOWN' buttons to scroll through the Program List. Upon selection you are returned to the Step Definition screen.

INS
STEP

LINE
↑

LINE
↓

PAGE
↑

PAGE
↓

Stop

The Stop Command is used to instruct the pump to stop operation. This Command can be used in conjunction with I/O Events to control pump operation based on external events.

TO INSERT A STOP COMMAND IN A PROGRAM:

1. From the Program Main screen, choose Step Definition to enter the Step Definition screen.
2. Highlight (select once) the Program Step that you want the Stop Command to follow.
3. Choose the 'INS STEP' button to access the Program Steps Menu screen.
4. Choose Advanced Options to enter the Set Advanced Options screen.
5. Choose Stop to insert the Stop Command into the Program.

INS
STEP

I/O EVENT TRIGGERS

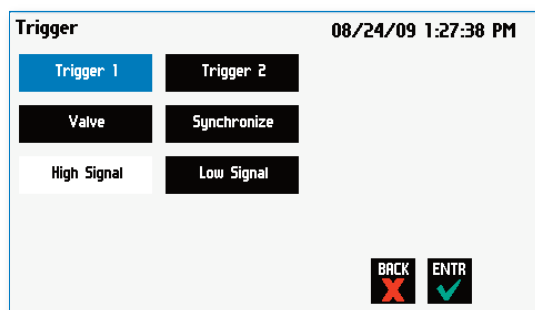
The Cole-Parmer® Touch Screen Series pump permits easy connection to external I/O devices via a variety of interfaces, including a built-in 15-pin D-sub connector. The Cole-Parmer® Touch Screen Series software incorporates simple commands for incorporating I/O controls into your pump Programs. Integrated commands for control of an external valve, multi-pump synchronization, and general trigger signals are integrated in the software. The following procedure describes the use of the Advanced Options I/O commands. Following the procedure is a description of the selection menu options for each of the I/O command types. Details for the electrical connection controls and pin-out description for the 15-pin User I/O Connector are provided in Appendix D of this manual.

The Cole-Parmer® Touch Screen Series manages trigger events through specific pins on the I/O Connector port on the back of the pump. The following details define pin logic:

- Triggers are communicated via Pin 4 or Pin 5.
- Valve control signals are communicated via Pin 8.
- Pump synchronization commands are communicated via Pin 6.
- Event signals are communicated via Pin 2.

I/O (OUT) Trigger

The I/O (OUT) Trigger command is used to send a signal to an external device when the Step in the Program is reached.



I/O (OUT) Trigger Advanced Options menu

The following options are available when setting an I/O (OUT) Trigger command:

- **Trigger 1:** Specifies output via the Trigger 1 (Pin 4) port. When selected, you must also specify whether the signal to be sent is High Signal or Low Signal.
 - **High Signal:** Sets the output signal to logic high.
 - **Low Signal:** Sets the output signal to logic low.
- **Trigger 2:** Specifies output via the Trigger 2 (Pin 5) port. When selected, you must also specify whether the signal to be sent is High Signal or Low Signal.
 - **High Signal:** Sets the output signal to logic high.
 - **Low Signal:** Sets the output signal to logic low.
- **Valve:** Used to trigger an external flow control valve via Pin 8. When selected, you must also specify the trigger action: Off or Valve Actuated.
 - **Valve Off:** Outputs a logic low to set valve to the off position.
 - **Valve Actuated:** Outputs a logic high to set valve to the actuated (open) position.
- **Synchronize:** Select this to send a signal via Pin 6 to a connected pump to trigger activity. When selected, you must define whether the desired action is to start an Infuse activity or a Refill (Withdraw) activity.
 - **Start Infuse:** Defines the rising edge signal to trigger an infusion action.
 - **Start Withdraw:** Defines the falling edge signal to trigger a refill (withdraw) action.

TO INSERT AN I/O (OUT) TRIGGER COMMAND:

1. From the Program Main screen, choose Step Definition to enter the Step Definition screen.
2. Highlight (select once) the Program Step that you want the I/O (OUT) Trigger to follow.
3. Choose the **'INS STEP'** button to access the Program Steps Menu screen.
4. Choose **Advanced Options** to enter the Select Advanced Options screen.
5. Choose the **I/O (OUT) Trigger** button to enter the Trigger Step Definition Screen.
6. Set your Trigger settings by choosing the appropriate option from the following sub-procedures:
 - a. To send a Trigger signal to an external device via Pin 4, select the Trigger 1 button. To send a Trigger signal via Pin 5, select the Trigger 2 button (selection will highlight in blue). Then select whether the signal will be High or Low by selecting the appropriate button (selection will highlight in white).
 - b. To send a trigger signal to an external valve, select Valve. Then select whether the trigger will actuate (Actuated button) or close (Off button) by selecting the appropriate button.
 - c. To send a synchronization signal to a connected pump, select Synchronize. Then select whether you want to start an infusion action (Start Infuse button) or a withdraw action (Start Withdraw button) by selecting the appropriate button.
7. Choose the Accept button to save your changes and return to the Step Definition screen.

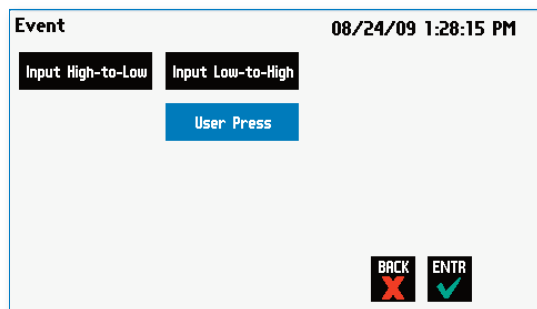
INS
STEP



NOTE: When using triggers, the first step in your program should define the trigger starting level (High or Low). You may then use the trigger command to change the level as desired for your Program.

I/O (IN) Event

The I/O (IN) Event command is used to have the pump pause and wait for a signal from an external device or for a manual user input at the touchscreen. These external signals are communicated via Pin 2 on the 15-pin I/O Connector port.



I/O (IN) Event Advanced Options menu

The following options are available when setting an I/O (IN) Event command.

- **Input High-to-Low:** Sets the pump to wait for a logic low input signal falling edge.
- **Input Low-to-High:** Sets the pump to wait for a logic high input signal rising edge.
- **User Press:** Sets the pump to wait for manual input from the user at the touchscreen (screen touch).

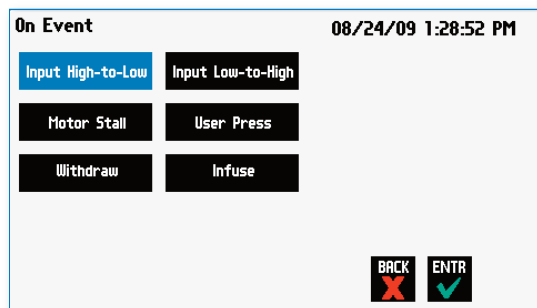
TO INSERT AN I/O (IN) EVENT COMMAND:

1. From the Program Main screen, choose Step Definition to enter the Step Definition screen.
2. Highlight (select once) the Program Step that you want the I/O (IN) Event to follow.
3. Choose the **'INS STEP'** button to access the Program Steps Menu screen.
4. Choose **Advanced Options** to enter the Select Advanced Options screen.
5. Choose the **I/O (IN) Event** button to enter the Event Step Definition Screen.
6. Select either Input High, Input Low or User Press, as appropriate.
7. Choose the Accept button to save your changes and return to the Step Definition screen.

INS
STEP

I/O (IN) Event Link

The I/O (IN) Event Link command is used to have the pump jump to a linked program if the signal is received. The external signals are communicated via Pin 2 on the 15-pin I/O Connector.



Advanced Options I/O (IN) Event Link Selection Screen

The following options are available when setting an I/O (IN) Event Link command. This command requires that you have a program set with the specific actions (steps) you want to perform if the event occurs so that you can choose the Program name being linked by this command.

- **Input High-to-Low:** Will link to the desired program upon a logic low input signal falling edge.
- **Input Low-to-High:** Will link to the desired program upon a logic high input signal rising edge.
- **User Press:** Will link to the desired program upon a manual input.
- **Motor Stall:** Will link to the desired program upon a Motor Stall event.
- **Withdraw:** Will link to the desired program upon a Withdraw activity.
- **Infuse:** Will link to the desired program upon an Infuse activity.

TO INSERT AN I/O (IN) EVENT LINK COMMAND:

1. From the Program Main screen, choose Step Definition to enter the Step Definition screen.
2. Highlight (select once) the Program Step that you want the I/O (IN) Event Link command to follow.
3. Choose the **'INS STEP'** button to access the Program Steps Menu screen.
4. Choose **Advanced Options** to enter the Select Advanced Options screen.
5. Choose the **I/O (IN) Event Link** button to enter the On Event Step Definition Screen.
6. Select the desired Event type that you want to have the event associated with from the options available (Input Low, Input High, User Press, Motor Stall, Withdraw, Infuse).
7. After selecting the desired option in the On Event screen, choose Accept to access the Program List selection screen. Choose the Program name to jump to. To use this command, you must have defined a linkable Program.
8. Choose the Accept button to save your changes and return to the Step Definition screen.

**INS
STEP**



NOTE: Since the 'on event' will cause a jump (link) to a different program, be cautious not to use embedded links in your program programming.

MANAGING COLE-PARMER® TOUCH SCREEN SERIES PROGRAMS

All Existing and User-Defined Programs are stored in the Program List for ongoing use. Over time, you may want to manage previously created Programs to keep your Program List easy to navigate. In addition, the Cole-Parmer® Touch Screen Series allows you to export Program information to an externally connected pump or to a computer that will be used to control operation. This section of the Guide describes the file management facilities available on the Cole-Parmer® Touch Screen Series.



NOTE: The program management procedures outlined in this section are not available for basic pump modes (Infuse Only, Withdraw/Infuse, etc.).

Saving A Copy Of A Program

Follow these instructions to save a copy of a Program under a new name:

1. From Program Configuration screen, choose the Program Name button to enter the Mode Selection screen.
2. From the list of available Programs, select (press once) the Program you want to save a copy of, then choose the File Options button to enter the File Options screen.
3. Choose **Save As/Copy** to access the Program Name Entry screen. The current Program Name will be displayed in the text entry box.
4. On the Program Name screen, type the new name for your Program using the onscreen keyboard. Choose the Symbols/Numbers button (.?123) to display the list of available numbers and symbols that can be included in the Program Name. To switch the keyboard back to letters mode, choose the Letters button (ABC). [15 characters maximum]
5. When finished entering the new Program Name, choose Accept to save your changes and return to the File Options screen.
6. Choose **'BACK'** to return to the Mode Selection screen.



Renaming A Program

Follow these instructions to rename a previously saved Program:

1. From the Program Configuration screen, choose the Program Name button to enter the Program Mode screen.
2. From the list of available Programs, select (press once) the Program you want to rename, then choose the File Options button to enter the File Options screen.
3. Choose **Rename Program** to access the Program Name Entry screen. The current Program Name will be displayed in the text entry box.
4. On the Program Name screen, type the new name for your Program using the onscreen keyboard. Choose the Symbols/Numbers button (.?123) to display the list of available numbers and symbols that can be included in the Program Name. To switch the keyboard back to letters mode, choose the Letters button (ABC). [15 characters maximum]
5. When finished entering the new Program Name, choose Accept to save your changes and return to the File Options screen.
6. Choose **'BACK'** to return to the Program Selection screen.



Appending A Program

Follow these instructions to Append a previously created Program into a new Program:

1. From the or Program Configuration screen, choose the Program Name button to enter the Mode Selection screen.
2. From the list of available Programs, select (press once) the Program you want as the base Program.
3. Choose the **File Options** button to enter the File Options screen.
4. Choose **Append** to access the Program to Append screen.
5. Choose (press twice) the Program you would like to Append (attach) to your base Program. The pump will return you to the File Options screen and display a confirmation that the Program was copied.
6. Choose **BACK** to return to the Mode Selection screen.



Deleting A Program

Follow these instructions to delete a previously saved Program:

1. From the Program Configuration screen, choose the Program Name button to enter the Mode Selection screen.
2. From the list of available Programs, select (press once) the Program you want to delete, then choose the File Options button to enter the File Options screen.
3. Choose **Delete Program**. The application will display a warning message in the Message Area asking you to confirm your deletion request.
4. Choose **Delete Program** again to confirm the deletion and return to the Mode Selection screen.

COLE-PARMER® TOUCH SCREEN SERIES PROGRAM EXPORT AND IMPORT

The Cole-Parmer® Touch Screen pump permits the daisy chaining of multiple pumps, as well as integration with/control by an external computer through HyperTerminal. To facilitate operation in these modes, the Cole-Parmer® Touch Screen Series application software includes a variety of commands designed to simplify the export/import of Programs between the pump and external devices. This section outlines the use of these commands.



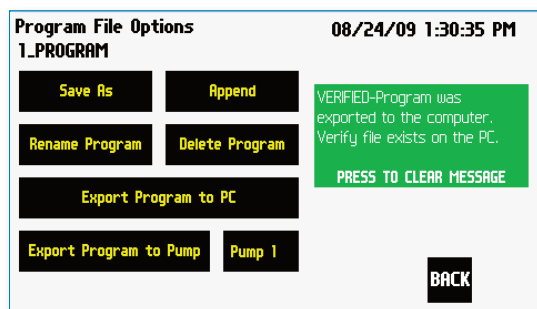
NOTE: Do not attempt to Export/Import Programs while the pump is running.



NOTE: Do not change screens while Importing/Exporting a method.



NOTE: For additional details on the physical interfaces used to connect the Cole-Parmer® Touch Screen Series to external devices, refer to the External I/O Connections section in the Appendices.



Export Program confirmation screen showing successful export

TRANSFER A PROGRAM BETWEEN PUMPS

Follow these instructions to transfer a Program between two Cole-Parmer® Touch Screen Series pumps.

1. Connect two pumps together via a RS-485 (firewire type) cable and set the address on the first pump to '00' using the Set Pump Address command in the Settings menu.
2. Set the address on the second pump to something other than '00' using the Set Pump Address command in the Settings menu.
3. From the Program Main screen, choose Program Select then choose the Program you would like to export.
4. Choose the **File Options** button.
5. Choose the Pump button and enter the address for the pump you would like to receive the Program.
6. Choose Export Program to Pump.
7. The pump will display a confirmation message in the Message Area on the screen when the export is complete. Note that long Programs may take a few seconds to export.
8. Before disconnecting the pumps, verify that the Program was successfully exported by choosing Program Select from the Program Main screen on the destination pump, then browse the available Programs to verify that the exported Program is shown in the menu.



NOTE: Before downloading or transferring a method to the pump, verify the number of steps available on the pump information screen.



Exporting A Program To An External Computer

Follow these instructions to export a Program to an external computer.

1. Connect the pump to the external computer via a RS-232 or USB cable.
2. On the computer, open the HyperTerminal application and verify that the computer can communicate with the pump.
3. At the pump, from the Program Main screen, choose Program Select then choose the Program you would like to export.
4. Choose the **'FILE OPTIONS'** button.
5. In HyperTerminal, choose Transfer→Capture Text...
6. Enter the File Name that you want to save the Program under, then choose Start.
7. Choose Export Program to PC.
8. The pump will display a confirmation message in the Message Area on the screen when the export is complete. Note that long Programs may take a few seconds to export.
9. At the destination computer, in HyperTerminal, when the Program is done transferring, the message "mend" will appear onscreen.
10. In HyperTerminal, choose Transfer→Capture Text...→Stop. Your Program is now saved on the external computer.

IMPORTING A PROGRAM FROM AN EXTERNAL COMPUTER

Follow these instructions to import a Program from an external computer.

1. Connect the pump to the external computer via a RS-232 or USB cable and set the Address value to '00' using the Set Pump Address command in the Settings menu.
2. On the computer, open the HyperTerminal application and verify that the computer can communicate with the pump.
3. In HyperTerminal, click on the Properties icons and choose the Settings tab on the pop-up box that is displayed.
4. Click on the ASCII Setup button.
5. In the Line Delay field, enter 500. This will set a 500 millisecond delay after each line is transferred during the import process.
6. Click OK twice to return to the HyperTerminal Main Menu screen.
7. In HyperTerminal, choose Transfer→Send Text File...
8. Select the Program that you want to import to the pump and click Open.
9. Before disconnecting, verify that the Program was transferred by choosing Program Select from the Program Main screen at the destination pump, then browse the Program menu to confirm that the new Program is now saved at the pump.

EXTERNAL PUMP CONTROL

This section of the Cole-Parmer® Touch Screen Series manual describes the control of the pump using an external computer device. If communicating via the USB Port, see Appendix D for one time installation instructions of USB Driver.

Setting Up HyperTerminal

HyperTerminal is a Windows application designed to support the external control of devices such as the Cole-Parmer® Touch Screen Series pump through an RS-232 or USB connection. The following instructions describe the configuration of the HyperTerminal application.

1. Select "Start - All Programs - Accessories - Communications - HyperTerminal".
2. Enter a name for a New Connection (i.e CP_Touchscreen), then click OK.

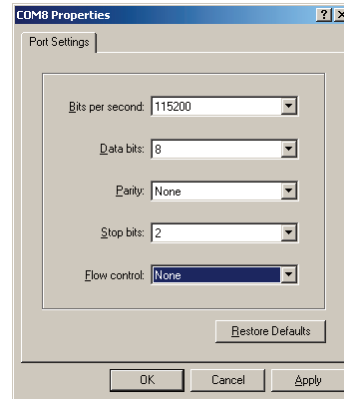


3. Select the Virtual Comm Port from the "Connect using" drop-down list. Click OK. (If the Virtual Comm Port is not known, use Device Manager to find it. Instructions on setting up the Virtual CommPort Driver are supplied in Appendix E.)




NOTE: To determine the proper comport, in your PC's "Systems Properties" screen, under the "Hardware" tab, select "Device Manger". Expand the "ports (COM & LPT)" section and you should see "Syringe Pump". The comport is listed next to this in parenthesis.


4. Set up the Port Settings as shown below and click OK.



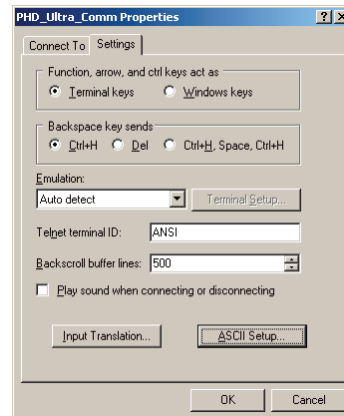5. Verify the Settings are as shown below.



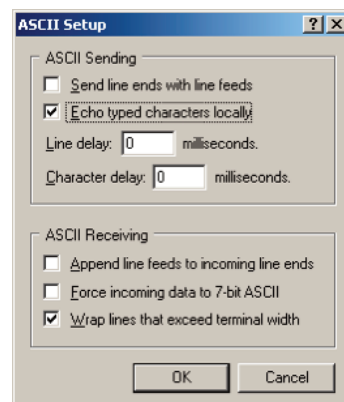6. Choose ASCII Setup and select "Echo typed characters locally" then click OK to complete the setup.



NOTE: To determine if your pump is connected correctly, in Hyperterminal's main window press 'Enter'; a ':' should appear. If not connected properly, close Hyperterminal, turn off the pump, unplug then re-plug the USB cable and re-open Hyperterminal.

PUMP CHAIN COMMANDS

The Pump Chain commands allow all pump control information to be managed from an external computer source. These commands can also be used to control a series of pumps (up to 100) from a single computer interface.

Cole-Parmer® Touch Screen Series commands are communicated to the pump via the RS-232 or USB port interfaces through a terminal program such as HyperTerminal. In using the Pump Chain commands, you will need to assign each pump in the pump chain a unique address, using the Set Pump Address command in the Settings menu. The address range is from 00 to 99. This address value is used to identify which pump is to receive a command and which pump is responding. The first pump in the chain, the one connected to the computer, must use address 0 (zero). Configure each pump with its assigned address and baud rate as described in the Pump Settings section of this guide.



NOTE: When daisy-chaining, all pumps must use the same baud rate. Values above 38400 are recommended.



NOTE: Once communication is established, if you manipulate the touchscreen, you must reestablish communications to the PC (<CR> will do this)



NOTE: System commands and start/stop commands can be executed from program or mode configuration screens. For parameter commands (rate, dia, vol, etc) it is necessary to execute from the mode configuration screen.

Using the Cole-Parmer® Touch Screen Series Pump Chain Commands

The following instructions will help you to utilize the Pump Chain commands feature on the Cole-Parmer® Touch Screen Series pump.

Commands may be abbreviated to the first four letters, i.e. address would be abbreviated addr. A space must follow the command if arguments are included.

If the pump address is nonzero, the one or two-digit pump address precedes the command. For example, to set the infuse rate for pump 12, the command would look like "12irat 3.2 ul/min".



NOTE: If the USB port is used on a pump with a nonzero address, the pump address does not need to be used.

Every line of the response from a pump with a nonzero address is prefixed by the two digit pump address (prefixed with a 0 if necessary) followed by a colon. No colon is used to separate the pump address from the prompt however.

In the command list below, the following convention is used:

- { } A parameter (without the brackets)
- [] Optional parameter
- | Separator between parameter choices
- A range of values

The following prompts are returned after a command is executed:

- : The pump is idle
- > The pump is infusing
- < The pump is withdrawing
- * The pump stalled
- T* The target was reached

If the pump is in poll mode, an XON character is added after the prompt.

System Commands

ADDRESS

Sets or displays the pump address. Valid range is 0 to 99.

address [0-99]

BAUD

Sets or displays the RS-232 port baud rate. Valid baud rates are 9600, 19200, 38400, 57600, 115200, 230400, 460800 and 921600. Note that some computers may not be able to handle baud rates above 115200.

baud [9600 | 19200 | 38400 | 57600 | 115200
230400 | 460800 | 921600]

DIM

Sets or displays the backlight level in percent. Valid range is 0 to 100. Note that specifying 0 will turn the backlight off so that the screen will not be visible.

dim [0-100]

ECHO

Sets or displays the RS-232/USB echo state. Valid states are on or off.

echo [on | off]

FREE

Displays the number of free program steps left on the disk.

free

FORCE

Sets or displays the infusion force level in percent. Valid range is 1 to 100.

force [1-100]

POLL

Sets or displays the polling mode state. In polling mode, prompts are inhibited if they aren't generated by a typed command. If polling mode is set on the master, an XON character is issued after every prompt. Valid states are on or off.

poll [on | off]

TIME

Sets or displays the date and time.

time [mm/dd/yy] [hh:mm:ss]

TILT

Calibrates the tilt sensor.

tilt



NOTE: Pump must be oriented in the horizontal plane during tilt calibration.

VALVE

Sets or displays the valve state. Valid states are on or off.

valve [on | off]

VDUTY

Sets or displays the valve duty cycle in percent. Valid range is 1 to 100.

vduty [1-100]

VER

Displays the short version string.

ver

VERSION

Displays the full version string.

version

Run Commands**IRUN**

Runs the pump in the infuse direction.

irun

WRUN

Runs the pump in the withdraw direction.

wrun

RRUN

Runs the pump in the opposite direction.

rrun

RUN

Simulates a key press of the run button located in the lower right corner of the screen.

run

STOP / STP

Stops the pump.

stop

stp

Rate Commands**CRATE**

Displays the current rate that the motor is running at. If the motor is stopped, the current rate is 0.

crate

DIAMETER

Sets or displays the syringe diameter in mm.

diameter [{syringe diameter}]

IRAMP †

Sets or displays the target infusion rate while ramping.

iramp [{start rate} {start units} {end rate}
{end units} {ramp time in seconds}] *

IRATE

Sets or displays the infusion rate.

irate [{rate} {rate units}] *

WRAMP †

Sets or displays the target withdraw rate while ramping.

wramp [{start rate} {start units} {end rate}
{end units} {ramp time in seconds}] *

WRATE

Sets or displays the withdraw rate.

wrate [{rate} {rate units}] *

* m, u, n, p/h, m, s

ex: m/m = milliliter/minute

n/s = nanoliter/second

† IRAMP & WRAMP commands do not show as changes to the mode displayed on the screen.

Volume Commands**CIVOLUME**

Clears the infused volume.

civolume

CTVOLUME

Clears the target volume.

ctvolume

CVOLUME

Clears both the infused and withdrawn volumes.

cvolume

CWVOLUME

Clears the withdrawn volume.

cwvolume

IVOLUME

Displays the infused volume.

ivolume

TVOLUME

Sets or displays the target volume.

tvolume [{target volume} {volume units}]
(units = m,u,n, p)

WVOLUME

Displays the withdrawn volume.

wvolume

Time Commands**CITIME**

Clears the infused time.

citime

CTIME

Clears both the infused and withdrawn times.

ctime

CTIME

Clears the target time.

ctime

CWTIME

Clears the withdrawn time.

cwtime

ITIME

Displays the infused time.

itime

TIME

Sets or displays the target time.

time [{target time} {time units}]
(units = h,m,s)

WTIME

Displays the withdrawn time.

wtime

Digital I/O Commands**INPUT**

Reads and displays the trigger input port status.

input

OUTPUT

Sets the level on one of the output ports.

output {1 | 2} {high | low}

SYNC

Sets the level on the sync port.

sync {high | low}

ERROR CODES

NOTE: Error and Alarm text are shown in a yellow box whereas Notice, Confirm and Verified are shown in a green box.

Run Mode Errors: Error 1XX

Set-up Errors: Error 2XX

Program Errors: Error 3XX

Critical Errors: Error 4XX

Run Mode Errors

ERROR 101- Power to the pump was lost.

ERROR 102- The A/D convertor is not working.

Set-up Errors

ERROR 201- Display brightness percentage is out of range, enter a value from 10 to 100

ERROR 202- Near end of run alarm percentage is out of range, enter a value from 1 to 100

ERROR 203- Force percentage is out of range, enter a value from 1 to 100.

ERROR 204- Pump address is out of range, enter a value from 0 to 99.

ERROR 205- Pump to pump address is out of range, enter a value from 1 to 99.

ERROR 206- Custom syringe volume is out of range, enter a value from 0.05 µl to 250 ml.

ERROR 207- Custom syringe diameter is out of range, enter a value from 0.1 mm to 99 mm.

ERROR 208- Date is out of range, enter a month from 1 to 12, a day from 1 to 31 and a year from 08 to 37.

ERROR 209- Time is out of range, enter the hour from 1 to 12 (or 0 to 23) and min and sec from 0 to 59.

ERROR 210- Infuse Rate is out of range.

ERROR 211- Withdraw rate is out of range.

ERROR 212- Syringe diameter is out of range.

ERROR 213- Target volume is out of range.

ERROR 214- Target time is out of range.

ERROR 215- Total volume exceeds syringe capacity.

ERROR 216- Volume is less than the minimum value.

ERROR 217- Motor has stalled. use the << or >> Key to move the pusher block. Check the force setting.

ERROR 218- Program needs a rate, enter a rate.

ERROR 219- Program needs infuse and withdraw rates and a target volume, enter these parameters.

ERROR 220- Date/Time is out of range.

Program Errors

ERROR 301- Delay time out of range, minimum value is 0.20 seconds.

ERROR 302- Number of steps must be between 2nd 500.

ERROR 305- Remote pump is busy or encountered and error. Retry and check remote pump.

ERROR 306- Download to remote pump not complete. Check remote pump.

ERROR 307- Syringe size and flow rate not compatible, cannot set motor speed.

ERROR 308- Program name already exists. Rename your program.

ERROR 309- Linked to program does not exist.

ERROR 310- Not all fields have been set. Check the configuration.

ERROR 311- Rate is out of range.

ERROR 312- Time is out of range.

ERROR 313- Volume is out of range.

ERROR 314- First rate is out of range.

ERROR 315- Second rate is out of range.

ERROR 316- First time is out of range.

ERROR 317- Second time is out of range.

ERROR 318- First volume is out of range.

ERROR 319- Second volume is out of range.

ERROR 320- The # of pulses is 0, enter a number larger than 0.

ERROR 321- Bad memory sector in step #%. Step was deleted.

ERROR 322- Target volume exceeds syringe capacity in step #%. Change the target volume.

ERROR 323- Incompatible rates vs. syringe size in step #%. Change either the rate or the syringe size.

ERROR 324- The # of repeats is 0, enter a number larger than 0.

ERROR 325- "Go To" cannot link to itself. Select another step.

Critical Errors

ERROR 401-	Infuse rate could not be saved to disk. Contact factory.
ERROR 402-	Withdraw rate could not be saved to disk. Contact factory.
ERROR 403-	Program step could not be saved to disk. Contact factory.
ERROR 404-	Program could not be read from disk. Contact factory.
ERROR 405-	Program step could not be read from disk. Contact factory.

Notice

NOTICE-	The display is locked. Press the lock button twice to activate the touch screen.
NOTICE-	Deleting program.
NOTICE-	Importing new program: %s...
NOTICE-	Exporting step #%d...
NOTICE-	Checking program.
NOTICE-	Pressing the run button will start the program at the following date and time: %s.
NOTICE-	Program is on a time delay. To cancel the delay and start the pump, press here.
NOTICE-	Waiting for screen press.
NOTICE-	Waiting for screen press or event.
NOTICE-	Waiting for an event.
NOTICE-	Copying a program.
NOTICE-	Verify rates and target.
NOTICE-	Withdrawal complete, press to silence alarm.
NOTICE-	Infusion complete, press to silence alarm.

Confirm

CONFIRM-	Press the lock button again, to confirm you want to unlock the touch screen.
CONFIRM-	Press button again to delete this program.
CONFIRM-	Press button again to reset calibration date.
CONFIRM-	Press button again to reset lubrication interval.
CONFIRM-	Press button again to delete step.

Verified

VERIFIED-	Program was copied.
VERIFIED-	Program was exported to the computer. Verify file exists on the PC.
VERIFIED-	Added step #%d...
VERIFIED-	Program was successfully imported.

Alarm

ALARM-	Pump is due for calibration.
ALARM-	Pump is due for lubrication.
ALARM-	Pump is due for calibration and lubrication.

Appendices

APPENDIX A: SYRINGE VOLUME/DIAMETER REFERENCE TABLE

KD Scientific Stainless Steel

Size	Diameter
2.5 ml	4.851 mm
8	9.525
20	19.13
50	28.6
100	34.9

Air-Tite HSW Norm-Ject

Size	Diameter
1 ml	4.69 mm
2.5	9.65
5	12.45
10	15.9
20	20.05
30	22.9
50	29.2

Becton Dickinson 'Plasti-pak'

Size	Diameter
1 ml	4.699 mm
3	8.585
5	11.989
10	14.427
20	19.05
30	21.59
50	26.594
60	26.594

Terumo

Size	Diameter
1 ml	4.70 mm
3	8.95 mm
5	13
10	15.8
20	20.15
30	23.1
60	29.1

SGE Scientific Glass Engineering

Size	Diameter
5 µl	0.343 mm
10	0.485
25	0.728
50	1.03
100	1.457
250	2.303
500	3.257
1 ml	4.606 mm
2.5	7.284
5	10.301
10	14.567
25	23
50	27.5
100	35

Hamilton Glass - All Types

Size	Diameter
0.5 µl	0.103 mm
1	0.146
2	0.206
5	0.343
10	0.485
25	0.729
50	1.03
100	1.457
250	2.304
500	3.256
1 ml	4.608 mm
1.25	5.151
2.5	7.285
5	10.3
10	14.567
25	23.033
50	32.573
100	32.573

Cadence Science, Inc. Formerly Popper & Sons MICRO-MATE Glass

Size	Diameter
0.25 ml	3.47 mm
0.5	3.62
1	4.82
2	8.91
3	8.91
5	11.71
10	14.65
20	19.56
30	22.7
50	28.02
100	35.7

Becton Dickinson Glass-All Types

Size	Diameter
0.5 ml	4.64 mm
1	4.64
2.5	8.66
5	11.86
10	14.34
20	19.13
30	22.7
50	28.6
100	34.9

KD Scientific Glass

Size	Diameter
1 ml	4.80 mm
2	6.45
5	12.60
10	15.50
20	20.40
30	22.90
50	27.45
100	35.90
150	34.20

Sherwood-Monoject Plastic

Size	Diameter
1 ml	4.674 mm
3	8.865
6	12.600
12	15.621
20	20.142
35	23.571
60	26.568
140	37.948

Hoshi

Size	Diameter
1 ml	6.50 mm
2	9.10
3	10.00
5	12.60
10	15.10
20	20.45
30	22.50
50	25.60
100	34.00

Natsume

Size	Diameter
0.25 ml	2.5 mm
0.50	2.8
1	4.0
2	6.2
3	7.4
5	9.2

Top

Size	Diameter
1 ml	6.40 mm
2.5	9.30
5	13.10
10	15.3
20	21.0
30	23.0
50	29.0

Nipro

Size	Diameter
1 ml short	6.6 mm
1 ml long	4.7
2.5 ml	9.0
5	13.0
10	15.8
20	20.1
30	23.2
50	29.1

Suggested Force Level Settings for Common Syringes

Syringe Type/Material	Capacity	Force Setting	Common Manufacturer(s)
Plastic Syringes	≤20ml	50%	BD plastic, Sherwood, Airtite, Terumo
	>20ml	100%	
Stainless Steel	All	100%	KDS, Cole-Parmer
Glass/Glass	≤20ml	20%	KDS Glass, Cadence (Popper)
	>20ml	30%	
Glass/Plastic	≤1000µl	30%	Hamilton, SGE
	≤5ml	50%	
	>5ml	100%	


APPENDIX B: MIN/MAX FLOW RATES

*Nominal Minimum/Maximum Flow Rates for Various Syringes.
(Actual Limits will vary depending on syringe manufacturer)*

Syringe Size	Syringe ID	Minimum Rate	Units	Maximum Rate	Units
0.5 μ l	0.103 mm	3.12	pl/min	1.589	μ l/min
1 μ l	0.146 mm	6.18	pl/min	3.180	μ l/min
2 μ l	0.206 mm	12.301	pl/min	6.358	μ l/min
5 μ l	0.343 mm	33.96	pl/min	17.630	μ l/min
10 μ l	0.485 mm	67.72	pl/min	35.249	μ l/min
25 μ l	0.729 mm	153.42	pl/min	79.640	μ l/min
50 μ l	1.03 mm	306.24	pl/min	158.984	μ l/min
100 μ l	1.457 mm	612.72	pl/min	318.126	μ l/min
250 μ l	2.304 mm	1.533	nl/min	795.51	μ l/min
500 μ l	3.256 mm	3.06	nl/min	1.588	ml/min
1000 μ l	4.608 mm	6.129	nl/min	3.181	ml/min
1 ml	4.699 mm	6.373	nl/min	3.308	ml/min
3 ml	8.585 mm	21.272	nl/min	11.044	ml/min
5 ml	11.989 mm	41.485	nl/min	21.539	ml/min
10 ml	14.427 mm	60.073	nl/min	31.190	ml/min
20 ml	19.05 mm	104.74	nl/min	54.383	ml/min
30 ml	21.59	134.533	nl/min	69.852	ml/min
50 ml	26.594 mm	204.122	nl/min	105.985	ml/min
100 ml	35.7 mm	367.839	nl/min	190.992	ml/min
140 ml	37.948 mm	415.623	nl/min	215.803	ml/min

APPENDIX C: RACK ACCESSORIES/OPTIONS

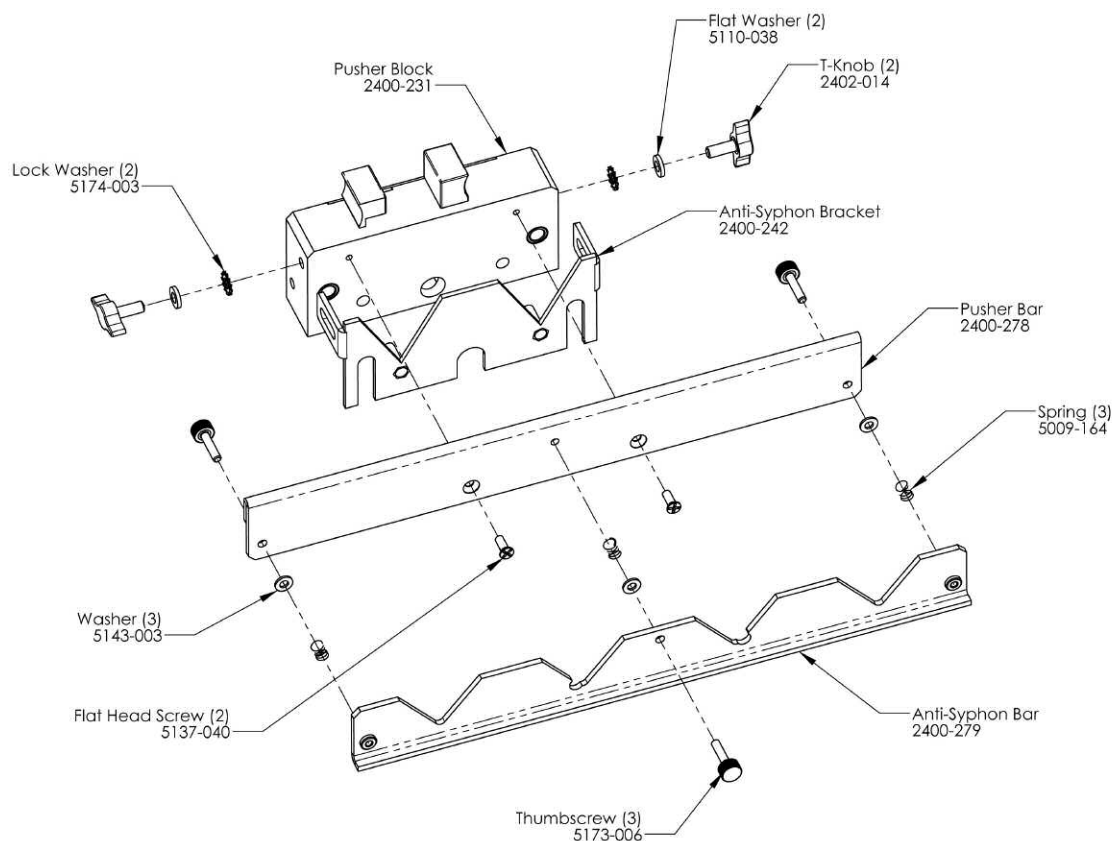
Mounting a 4 X 140 Multi-Rack (Continued on page 66)



Pusher Block Assembly


1. Remove Dual Syringe Anti-Syphon Bracket by unscrewing and removing the two T Knobs and two Spring Washers.
2. Attach the Pusher Bar to the Pusher Block using two flat head screws.
3. Attach the Anti-Syphon bar to the Pusher Bar with three Thumbscrews and three Springs. Slide Thumbscrews thru the Pusher Bar's holes, slide the Springs over the Thumbscrews and tighten Thumbscrews to the Pusher Bar.

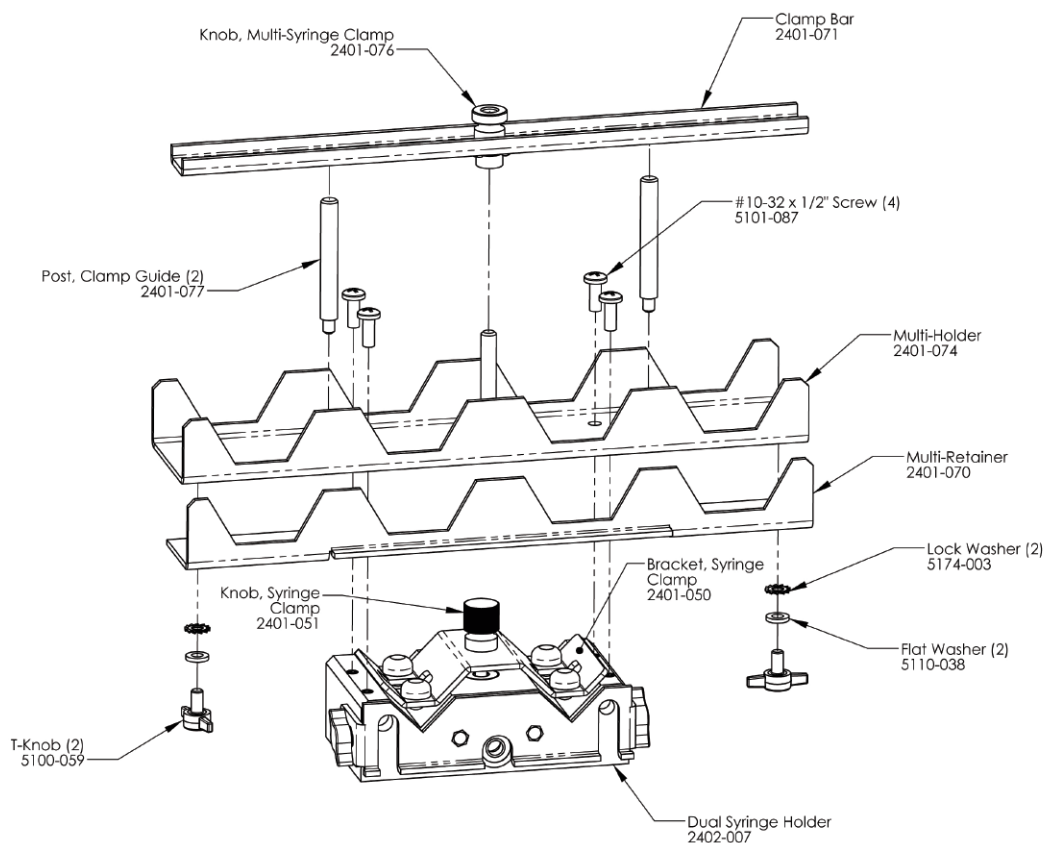Mounting a 4 X 140 Multi-Rack (Continued)




Syringe Holder Assembly

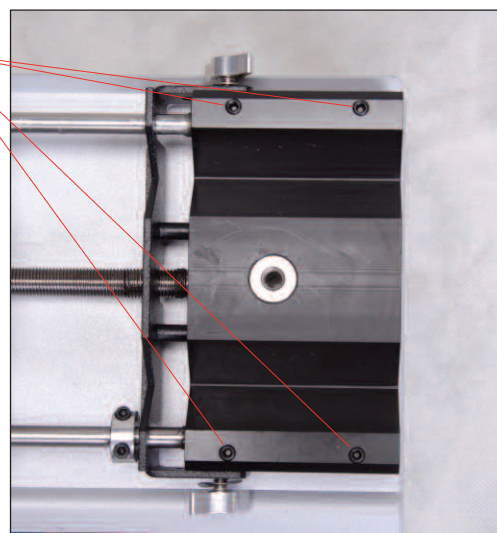1. Remove four set screws from the Dual Syringe Holder using the provided Hex Key.
2. Attach the Multi-Holder onto the Dual Syringe Holder using four #10-32 x 1/2" screws using a phillips head screwdriver.
3. Attach Clamp Bar by screwing the Multi-Syringe Clamp Knob into the threaded bar attached to the Multi-Holder.



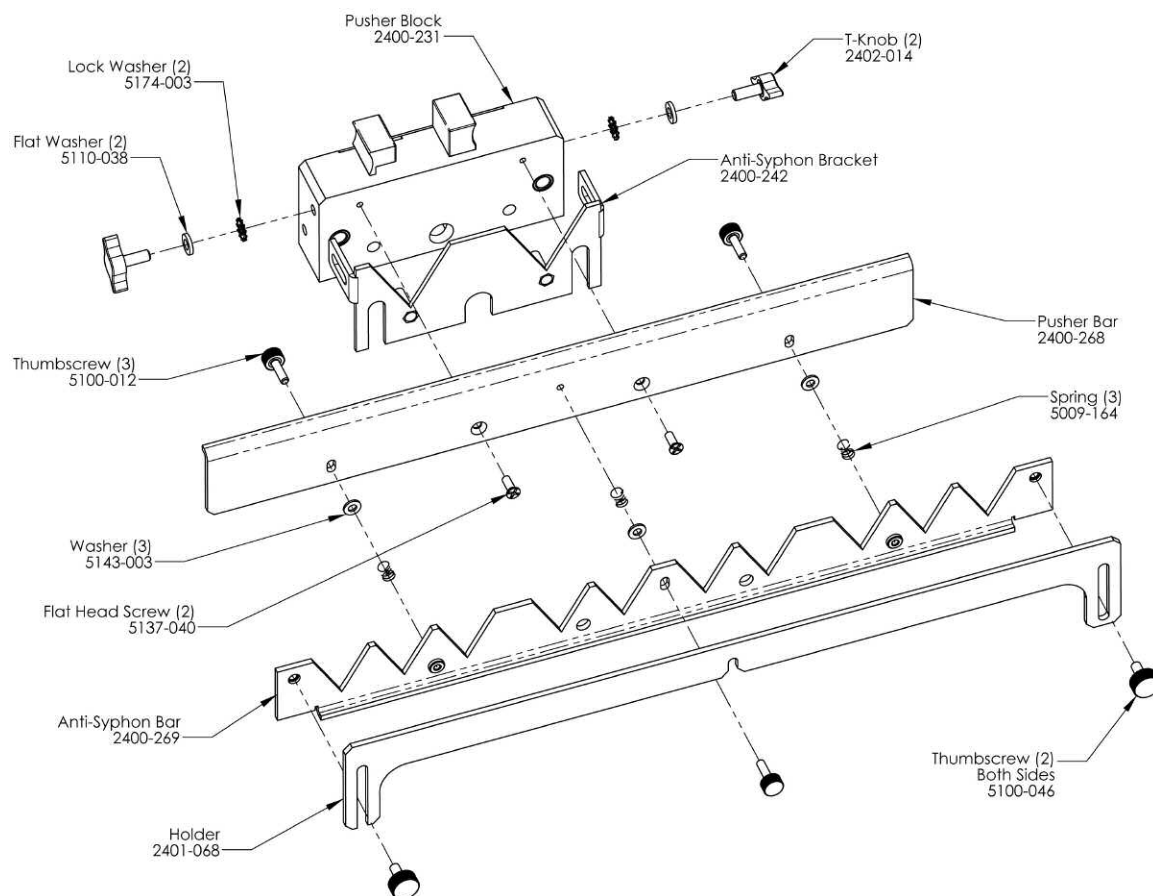
Mounting a 6/10 Multi-Rack (Continued on page 68)



Pusher Block Assembly

1. Remove Dual Syringe Anti-Syphon Bracket by unscrewing and removing the two T Knobs and two Spring Washers.
2. Attach the Pusher Bar to the Pusher Block with the two flat head screws using a phillips head screwdriver.

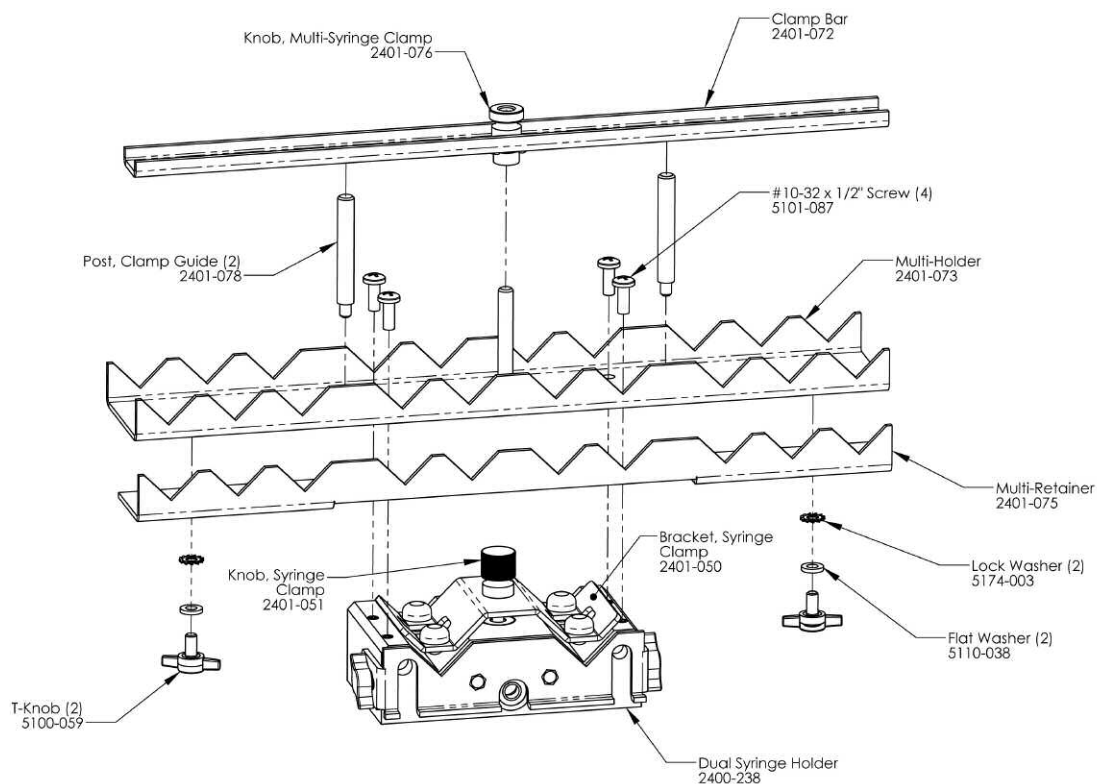Mounting a 6/10 Multi-Rack (Continued)




Syringe Holder Assembly

1. Remove Syringe Clamp bracket from the Syringe Holder by unscrewing knob.
2. Remove four set screws from the Dual Syringe holder using the provided Hex Key.
3. Place the Multi-Holder on the Dual Syringe Holder and attach using four #10-32 x 1/2" Screws using phillips head screw driver.
4. Screw the two Post Clamp Guides into Multi Syringe Holder
4. Attach Clamp Bar by screwing the Multi-Syringe Clamp Knob into the threaded bar attached to the Multi-Holder.

Syringe Loading - 4 x 140 & 6/10 Multi-Racks

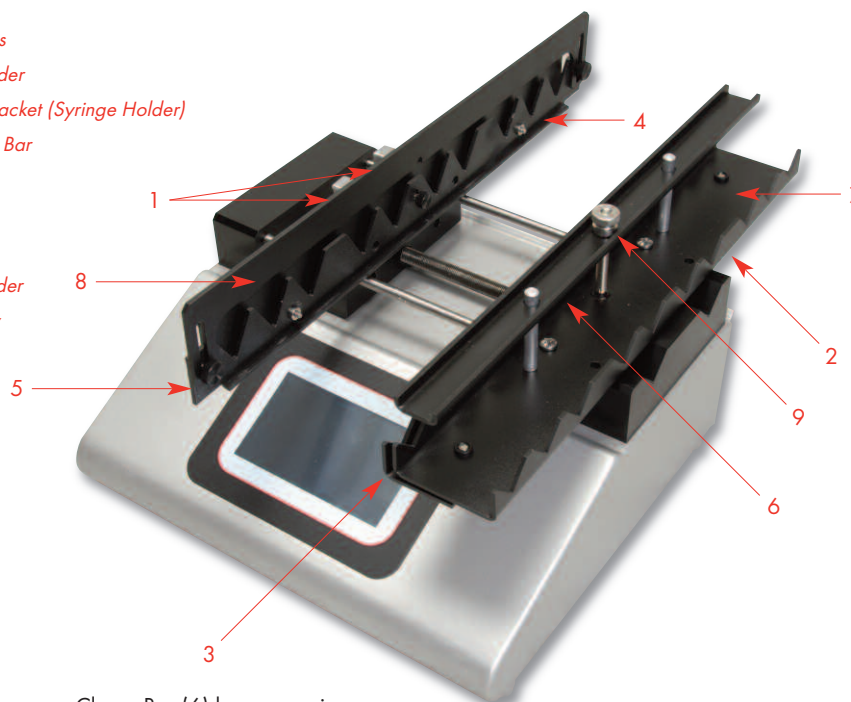1. Release Tabs
2. Syringe Holder
3. Retaining Bracket (Syringe Holder)
4. Anti-Syphon Bar
5. Pusher Bar
6. Clamp Bar
7. T-Knobs
8. Plunger Holder
9. Thumbscrew




1. Raise or remove Clamp Bar (6) by unscrewing thumb screw (9).
2. Loosen T-Knobs (7), located under the Syringe Holder (2), and slide the Retaining Bracket (3) on the Syringe Holder open.
3. Adjust the Pusher Block to the approximate length of the syringe(s) by squeezing the Release Tabs (1) and sliding the Pusher Block into position.
4. Loosen the three thumb screws on the Anti-Syphon Bar (4) to provide enough room for the plunger flange of the syringe(s). If using the 6/10 rack, loosen the two thumb screws on the outside of the Plunger Holder (8). Flip the Plunger Holder open.
5. Place the syringe(s) on the Syringe Holder (2) ensuring that the barrel flanges are within the Syringe Holder's Retaining Bracket (3) and the plunger flanges are between the Pusher Bar (5) and Anti-Syphon Bar (4).




NOTE: If using fewer syringes than the maximum number of syringes, try to place them symmetrically about the lead screw to present a balanced load.

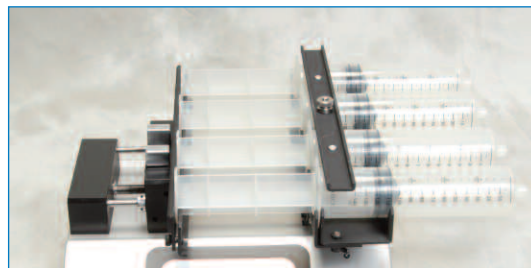6. Push the Syringe Holder's Retaining Bracket (3) tight against the barrel flanges and tighten into place using the two T-Knobs (7).
7. Tighten the Anti-Syphon Bar (4) using the three thumb screws. If using the 6/10 rack, flip the Plunger Holder back into place (8) and secure by tightening thumb screws.
8. Lower or replace the Clamp Bar (6) and tighten thumb screw until Clamp Bar is firm against syringe(s).




4 x 140 Multi-Rack

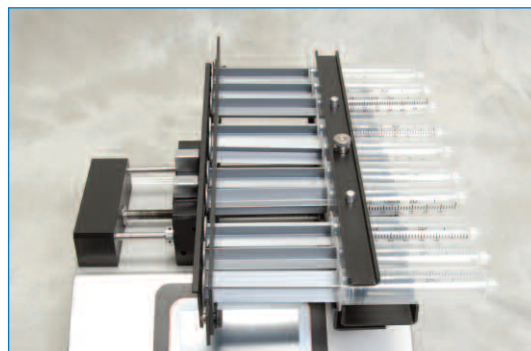

6/10 Multi-Rack

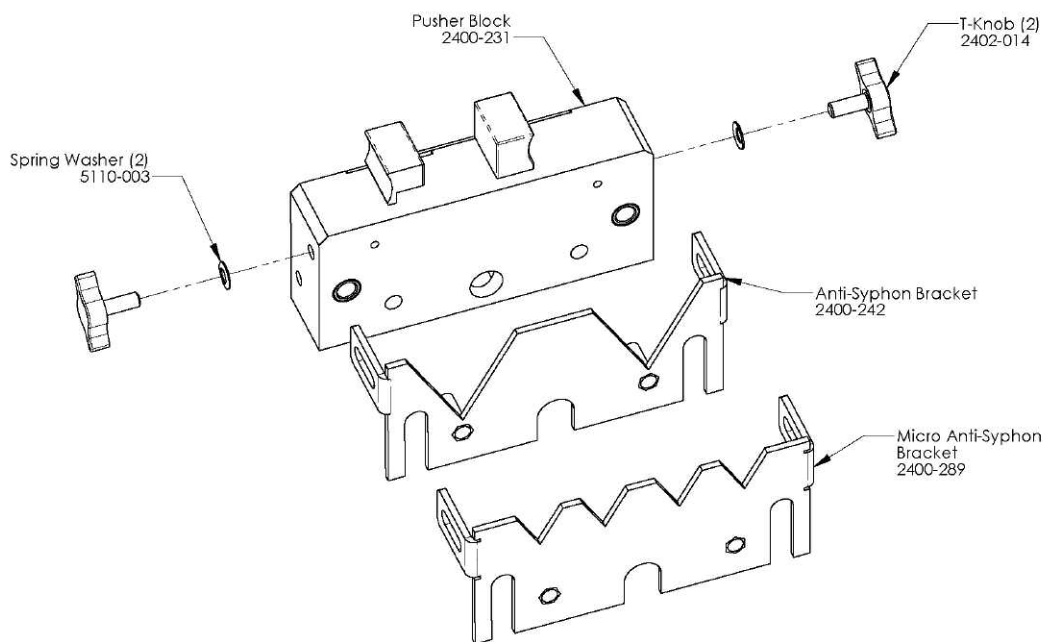
Mounting a Micro Liter Multi-Rack (Continued)



Pusher Block Assembly

1. Remove the Pusher Blocks Anti-Syphon Bracket by unscrewing and removing the T-Knobs (2) and Spring Washers (2).
2. Attach Micro Anti-Syphon Bracket using T-Knobs (2) and Spring Washers (2).

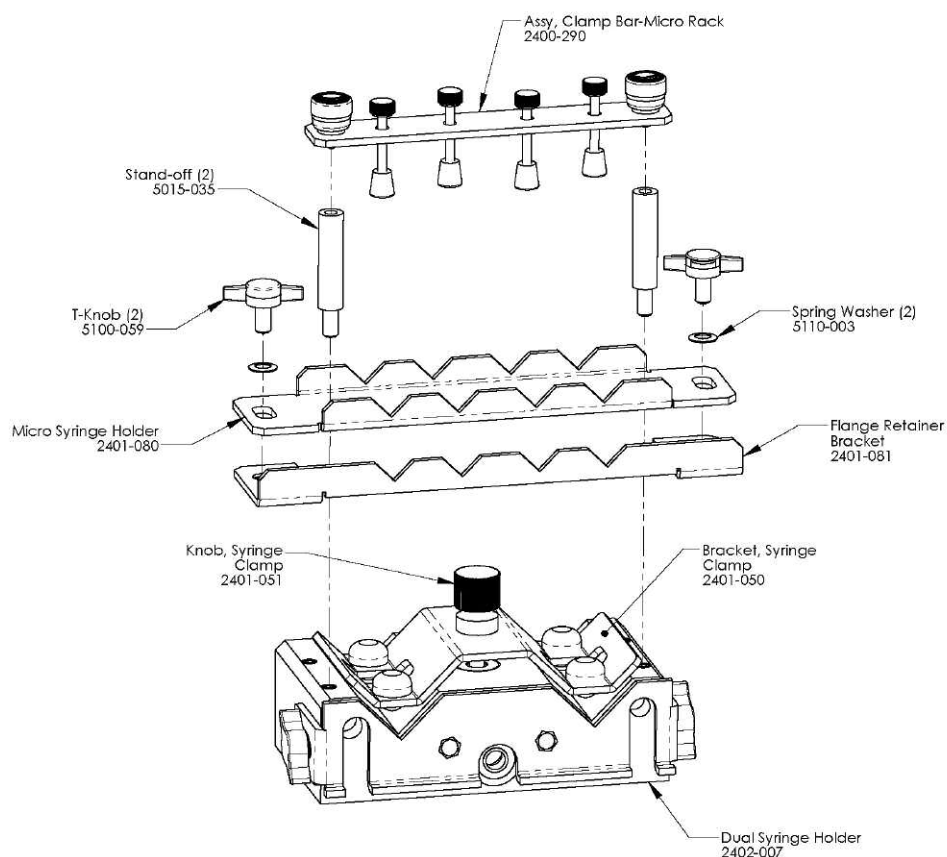Mounting a Micro liter Multi-Rack




Syringe Holder Assembly

1. Remove Syringe Clamp and knob from the Syringe Holder
2. Remove two set screws from the Dual Syringe Holder with the provided Hex Key.
3. Place the assembled Micro Syringe Holder onto the dual syringe holder and screw in the two Standoff Bars (2).
4. Place the Clamp Bar Assembly onto the Micro Syringe Holder and push and screw the Thumbscrew Caps (2) into the Standoff Bars.

APPENDIX C: RACK ACCESSORIES/OPTIONS

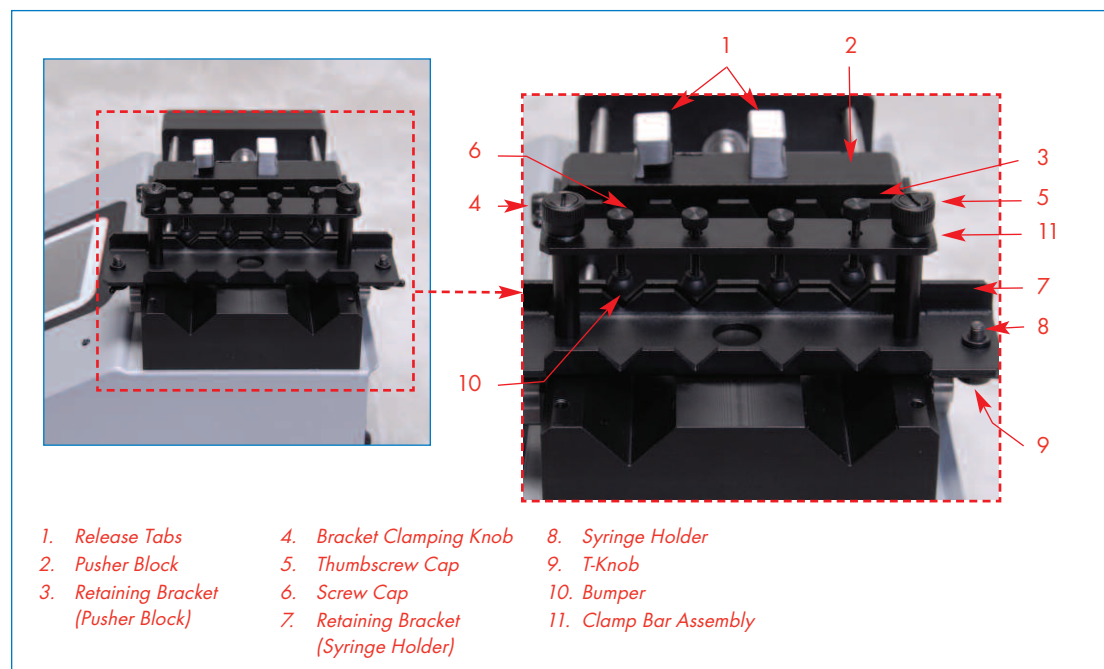Syringe Loading – Micro Liter Multi-Rack



Front view of Cole-Parmer® Touch Screen Series showing controls for micro liter multi-rack syringe loading

1. Squeeze Release Tabs (1) and adjust the Pusher Block (2) to the approximate length of the syringe(s).
2. Loosen the Retaining Bracket (3) on the Pusher Block (2) by unscrewing the Bracket Clamping Knobs (4).
3. Loosen the Retaining Bracket (7) on the Syringe Holder (8) by unscrewing the two T-Knobs (9).
4. Create room for the syringe(s) by unscrewing the Screw Caps (6) until the Bumpers (10) are at their maximum height.
5. Slide the syringe under the Bumper (10) placing the barrel flange within the Syringe Holder's Retaining Bracket (7) and the plunger flange within the Pusher Block's Retaining Bracket (3).
6. Lower the bumpers (10) until they are firm against the syringe(s) by screwing in the Screw Cap (6).
7. Push the Syringe Holder's Retaining Bracket (7) tight against the barrel flange and tighten into place using the T-knobs (9). Then push the Pusher Block's Retaining Bracket (3) tight against the plunger flange and tighten using the Bracket Clamping Knob (4).



NOTE: Mechanical stops are set using the same steps as a standard 2 syringe holder.

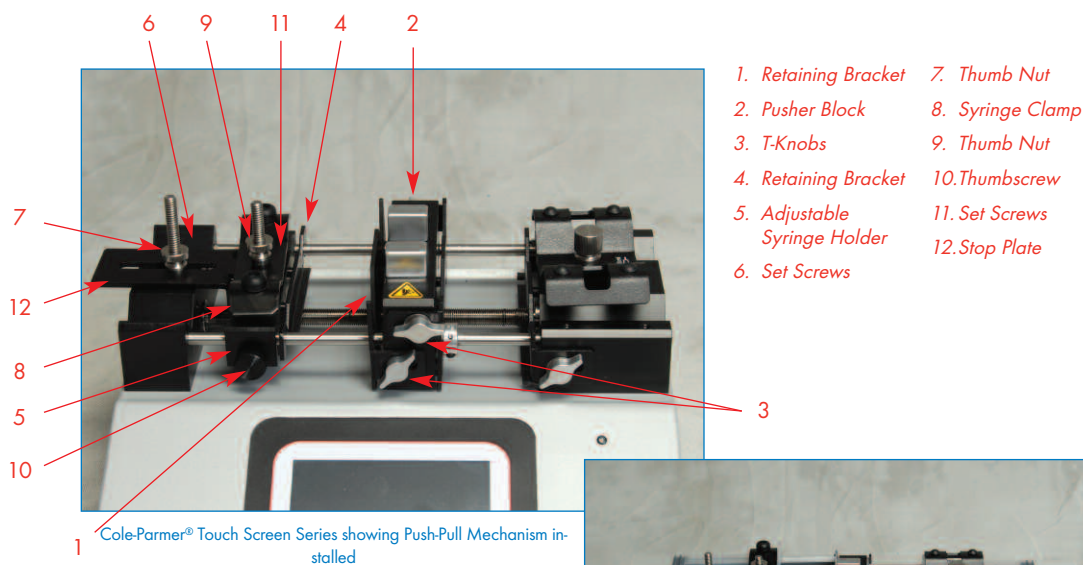


NOTE: If you are using syringes larger than 1 ml it may be easier to unscrew one of the Thumbscrew Caps (5) and swing the Clamp Bar Assembly (11) to the side, making more room to load the syringe(s). After the syringe(s) are in place swing the Clamp Bar back over the Syringe Holder (8) and reattach the Thumbscrew Cap (5).



Micro Liter Multi-Rack Configuration

Syringe Loading – Push-Pull Mechanism



Syringes < 60 ml

1. Loosen Thumb Nut (7) to release Adjustable Syringe Holder (5).
2. Loosen the Retaining Bracket (1) on the Pusher Block (2) using the two T-Knobs (3).
3. Remove Thumb Nut (9) and Syringe Clamp (8).
4. Loosen the Retaining Bracket (4) on the Adjustable Syringe Holder by unscrewing the two Set Screws (11) using an allen wrench.
5. Place the syringe(s) on the Adjustable Syringe Holder (5), extend as necessary then lock into place by tightening the Thumb Nut (7). Ensure that the syringe(s) flanges are within the Retaining Brackets.



6. Lock the syringe flanges into place by tightening the retaining brackets (1 & 4) on both the Pusher Block and Adjustable Syringe Holder using the Set Screws and T-Knobs.
7. Replace and secure the Syringe Clamp (8) and tighten using the Thumb Nut (9).



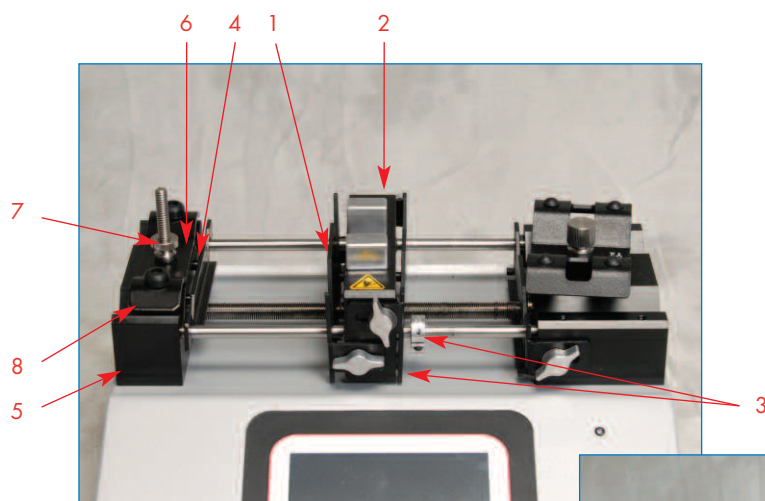
NOTE: For syringes <10ml you will need to remove the Thumb Nut (9) and then remove and flip over the Syringe Clamp (8).



NOTE: Push-Pull mechanism will support syringes of sizes 0.5µl to 140ml*

*Stroke length allows for only 120ml usable capacity from the 140ml syringe.

Syringe Loading – Push-Pull Mechanism (continued)



Syringes ≥ 60 ml

To remove Adjustable Syringe Holder:

1. Unscrew and remove both thumb nuts (7 & 9). Remove Syringe Clamp (8) and Stop Plate (12).
2. Loosen Set Screws (11) and remove Retaining Bracket (4).
3. Unscrew and remove thumbscrew (10) and then lift off the Adjustable Syringe Holder (5).
4. Attach the Retaining Bracket (4) to the remaining Syringe Holder using Set Screws (6).
5. Place the Syringe Clamp (8) on the Syringe Holder and replace Thumb Nut (7).



Loading Syringes ≥ 60 ml

1. Loosen Retaining Bracket (1) on the Pusher Block (2) using T-Knobs (3).
2. Loosen Retaining Bracket (4) on the Syringe Holder (5) using Set Screws (6).
3. Remove Thumb Nut (7) and Syringe Clamp (8). Place syringe(s) on the Syringe Holder (5), ensure that the syringe(s) flanges are within the Retaining Brackets.
4. Press the Retaining Bracket firmly against the barrel flange and tighten into place using Set Screws (6).
5. Replace Syringe Clamp (8) and tighten with Thumb Nut (7).
6. Push the Pusher Block's Retaining Bracket (1) firmly against the plunger flange and tighten using T-knobs (7).

APPENDIX D: EXTERNAL CONNECTIONS

User I/O Connector Specifications


Pin Assignments

1. **Direction control input**
Rising edge sets pump to infuse
Falling edge sets pump to refill
2. **Trigger Input**
Event trigger – falling edge triggers a program event
3. **Footswitch Input**
Momentary switch closure to ground or TTL Logic Low; falling edge toggles between states (run/stop)
4. **Trigger 1 Output**
TTL Logic Output-Default = Low Control thru program or serial comm
5. **Trigger 2 Output**
TTL Logic Output-Default = Low Control thru program or serial comm
6. **Sync Output**
TTL Logic Output –
Rising Edge = Start Infuse
Falling Edge = Start Refill
7. **Direction Output**
TTL Logic Output –
High = Infuse
Low = Refill
8. **Valve Output**
TTL Logic Output –
High = Valve Actuated
Low = Valve Off
- 9-13. **Signal Return / Ground**
14. **Run Indicator Voltage**
+5V through a 470 ohm resistor (connect to LED Anode)
15. **Run Indicator Output**
TTL Logic Output, active low (low = run) (connect to LED Cathode)

Electrical Specifications



Inputs:	$V_{IH} \geq 2V$	$I_{IH} \leq 20\mu A$
	$V_{IL} \leq 0.4V$	$I_{IL} \leq 0.5mA$
Outputs:	$V_{OH} \geq 3.8V$	$I_{OH} \leq 6mA$
	$V_{OL} \leq 0.4V$	$I_{OL} \leq 6mA$

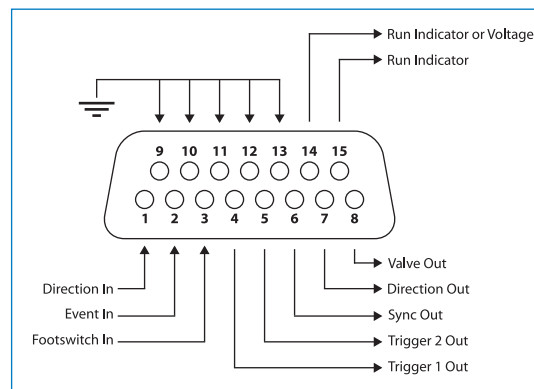All Inputs are pulled high to +5 Vdc through a 10k resistor All Outputs are at TTL Logic Levels




User I/O Connector Specifications

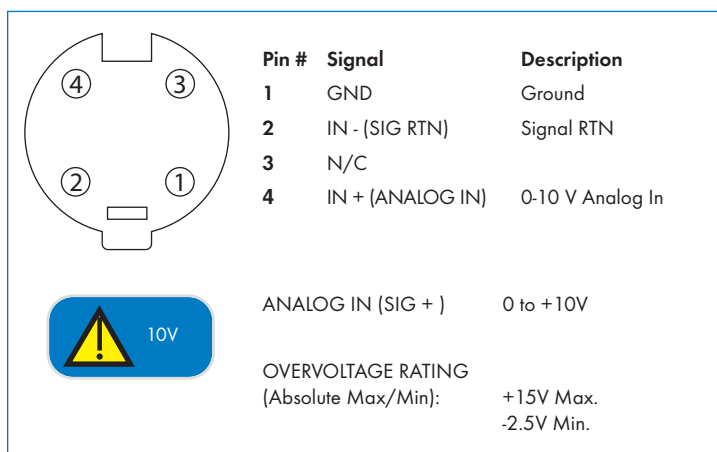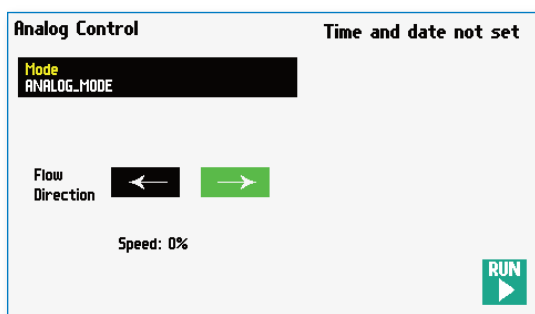Analog Input Option



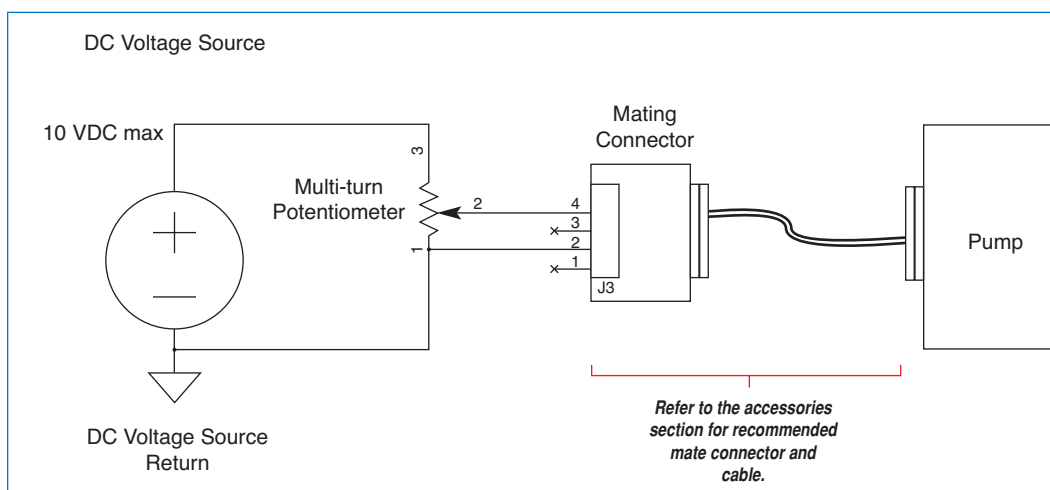Analog control of the motor speed can be achieved in many ways. The following circuit is recommended. A DC voltage (10 VDC max) will be needed to supply the excitation voltage to the circuit.




Connect circuit as shown in the above figure. Turning the Potentiometer CW will increase the motor speed (pin 4 is 10 VDC) and turning the Potentiometer CCW will decrease the motor speed (pin 4 is 0 VDC).



NOTE: Refer to the accessories section for recommended mate connector and cable.

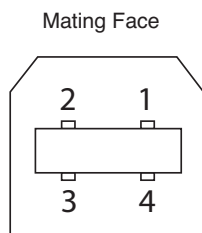Utilizing the USB Interface




Pin #	Signal
1	+5V
2	- Data
3	+ Data
4	GND

USB VIRTUAL COMMPORT DRIVER INSTALLATION

When you connect the Cole-Parmer® Touch Screen Series pump to a computer via USB for the first time, Windows will seek to install a driver for communication. The following section details the installation of the Virtual CommPort Driver supplied with the Cole-Parmer® Touch Screen Series Pump.




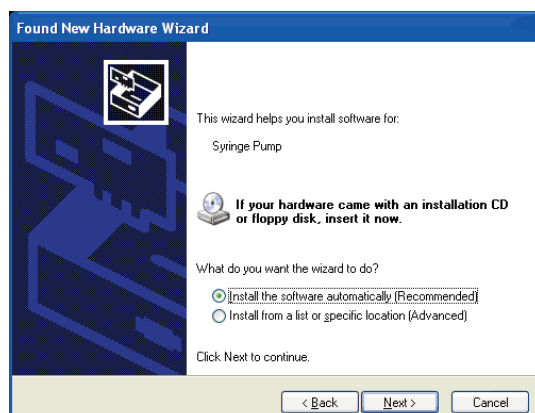1. On connection of the pump to the computer via USB, the Found New Hardware Wizard will launch and ask how you wish to locate the driver. Choose "No, not this time" and click Next.



2. Choose "Install from a list or specific location (Advanced)" and click Next.

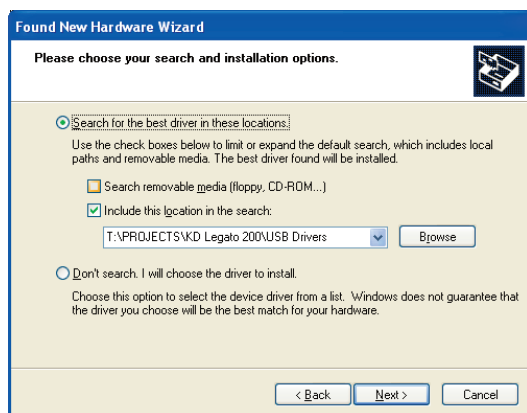
USB VIRTUAL COMMPORT DRIVER INSTALLATION (CONTINUED)




3. Browse for the proper USB driver installation folder in the USB folder and then click next. It can be found on the CD supplied with the Cole-Parmer® Touch Screen Series Pump.

USB Driver: WindowsXP or Vista

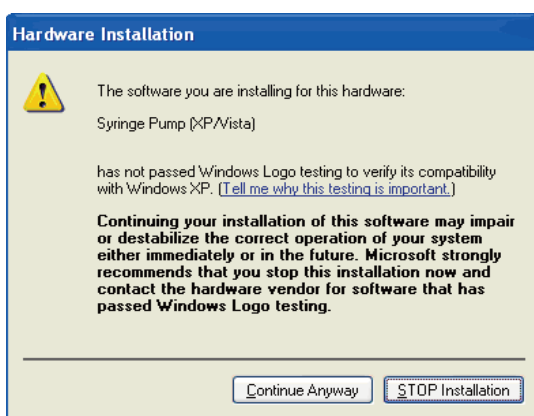USB Driver: Windows2000




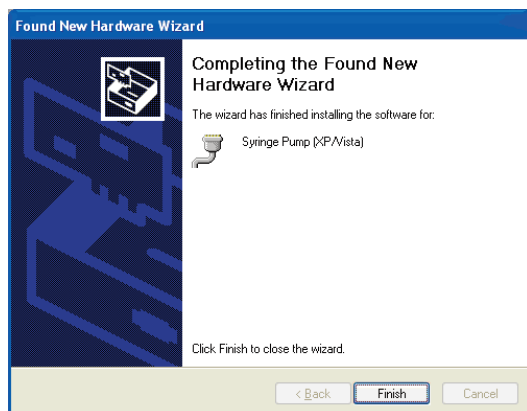4. If a warning message regarding Windows Logo testing is displayed, choose Continue Anyway to proceed with the installation.




5. The Virtual Commport software will be installed. When the software installation is complete, click Finish to close the wizard.

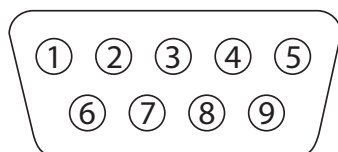RS-232 Specifications



Pin #	Signal
2	TXD
3	RXD
5	GND
1, 4, 6-9	N/C

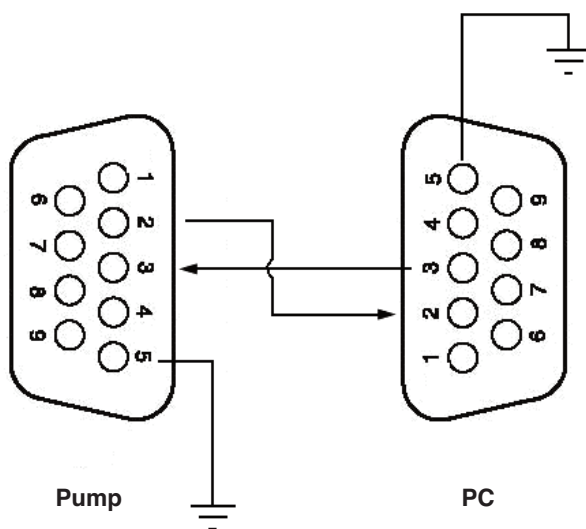RS-232 9-Pin D-SubMale



Baud Rates:

9600

19200

38400

57600

115200

Word size – 8

Parity – None

Stop Bits – 2

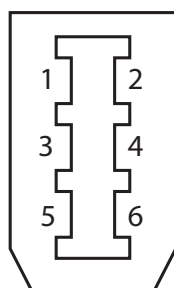
RS-232 Specifications

RS-485




CAUTION: DO NOT CONNECT TO FIREWIRE PORTS ON A PC. DAMAGE MAY OCCUR TO PUMP AND/OR PC.

Mating Face



Pin #	RS-485 In Signal	RS-485 Out Signal	Notes
1	PWR IN	PWR OUT (+30V)*	* Protected by resettable fuse
2	GND	GND	
3	RS-485 A	TIED TO PIN 4**	** 100 Ω , 1/2 w to GND termination
4	RS-485 B	TIED TO PIN 3**	
5	TIED TO PIN 6**	RS-485 A	
6	TIED TO PIN 5**	RS-485 B	

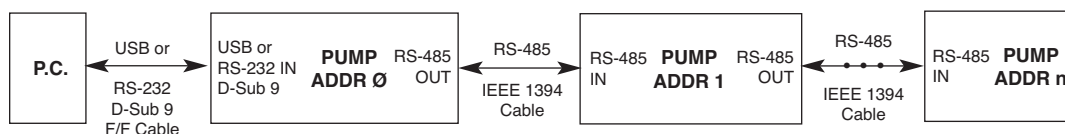RS-485 In/Out - IEEE-1394 Sockets



Daisy-Chaining via RS-485

APPENDIX E: MAINTENANCE

Maintenance

Cole-Parmer® Touch Screen Series pumps require no special maintenance other than keeping them clean by avoiding accidental spills of pumped material. Refer to the chemical compatibility chart below for specific chemical resistance information.

The two guide rods and the lead screw should be sparingly lubricated periodically; every 100 hrs the pump will remind you. The guide rods and the lead screw should be lubricated with Super Lube Synthetic Grease provided with the pump.

To clean the exterior surfaces, use a lint-free cloth to remove loose dust. Use care to avoid scratching the clear display window. For more efficient cleaning, use a soft cloth dampened [not soaked] with water, an aqueous solution of 75% isopropyl alcohol, or a mild detergent.


Chemical Compatibility Chart

Chemical	Concentration	Cover	Overlay Adhesive
Acetonitrile	99.5%	A	A
Methyl alcohol	99.8%	A	-
Hydrochloric Acid	36.5%	A	A
Bleach	100%	A	A
Windex®	100%	A	A
Isopropyl Alcohol	99.5%	A	D
Dish Soap	-	A	A
Motor Oil	100%	A	A
Simple Green	100%	A	A
Gasoline	100%	A	D
Ethylene Glycol (Antifreeze)	100%	A	A
Mineral Oil	100%	A	A
Ammonia	27.0%	A	A
Magnalube	100%	A	A
Superlube	100%	A	A
Acetone	99.5%	A	D
Soda (Coca-Cola)®	100%	A	A
Sulfuric Acid	95.0%	D	D
Dimethyl Sulfoxide	99.9%	A	C
Sodium Hydroxide	50%	A	A
Nitric Acid	68.0%	D	D
Chloroform	100%	D	D

- A Satisfactory
- B Fair, Minor Effects
- C Poor, Moderate Effects
- D Severe Effects
- not tested

The chemical resistance data provided is for reference only. The data is based on the manufacturer's recommendations and is not a guarantee of chemical compatibility for your test conditions. Wipe any spills immediately.

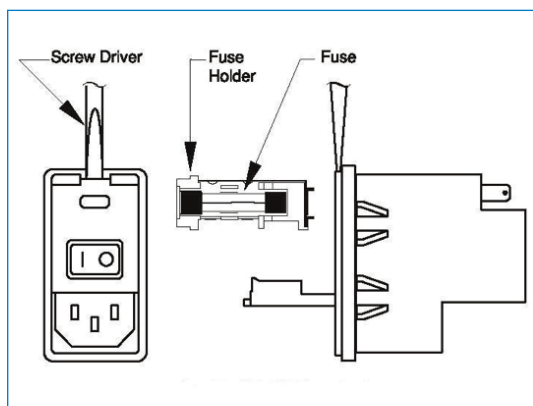Fuse Replacement



Key Information

1. Make sure the power cord is disconnected from the main supply before servicing the fuse.
2. Use only Type 3AG, 1/4" x 11/4, 1/2 amp, 250 volt; Type T (time delay) fuses (KD Scientific Part No. 5146-037 or equivalent).

Turn off power and remove power cord from power module. Use a straight blade screw-driver to pry open the access door. Remove the fuse holder and then remove the fuses from this holder as shown in figure below. Replace fuses, and then replace the fuse holder.

Battery Replacement



1. Remove screw and swing cover to the side.
2. Slide battery out from under clip.
3. Install Lithium Coin Battery: 3V, 16mm CR1620 (Part No. 5155-288 or equivalent) by sliding under clip with positive side facing out



CAUTION: USE ONLY SAME TYPE & RATED BATTERY. OBSERVE POLARITY WHEN INSTALLING.

Upgrading Cole-Parmer® Touch Screen Series Software

1. Upload the latest software version to your desktop (format is *filename.dfu*).
2. Connect the pump to PC using a USB cable.
3. On the set-up screen, press the Upgrade Software button twice to enter "Boot Mode". The pumps screen will be blank and the light should be flashing.

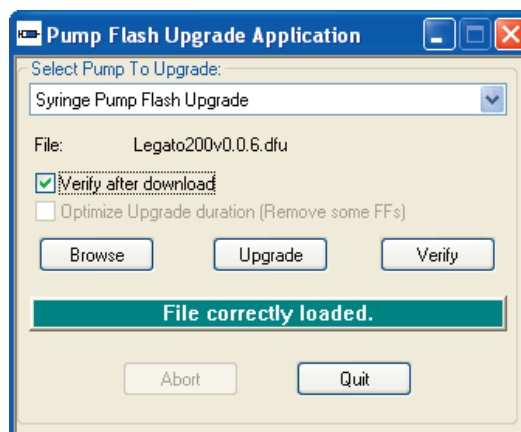


NOTE: To exit boot loader mode without upgrading the software power cycle the pump.

4. Disconnect then reconnect the USB cable.
5. From the CD provided with the pump, open the Cole-Parmer® Touch Screen Series Flash Upgrade Application.
6. Check "Verify after download" then click browse. Browse to the file previously saved on your desktop and click open.
7. Click upgrade. The previous software will first be erased, the new software will be loaded, and finally the software will be verified.
8. Disconnect the USB cable and power cycle the pump.



NOTE: First time users will need to install the Cole-Parmer® Touch Screen Series USB driver as well as the bootloader driver file (Cole-Parmer® Touch Screen Series-Bootloader.inf). Reference "USB Virtual Compact Driver Installation".



TROUBLESHOOTING

1. Infusion Accuracy

To ensure infusion accuracy always use new syringe(s) and measure syringe bore diameter and enter actual dimensions in millimeters (mm) using the Custom Syringe entry option. Additionally, make sure that the guide rods and lead screw are properly lubricated.

2. RS-232 Difficulties

Verify that the baud rates, data framing parameters, data bits, parity, and stop bits on all devices are the same. Verify flow control (handshaking) is set to None. If an address is set to other than address 0, make sure commands and queries are prefixed with the address. Make sure straightthrough pinning is used for the cables; do not use null-modem cables.



NOTE: It is recommended that you save your programs periodically by uploading them to a PC.

Cole-Parmer® Touch Screen Series Display Error Messages

3. Motor Stalled

This indicates that the pusher block travel has been impeded. This may be caused by the syringe plunger hitting bottom, a kink in the tubing (occlusion), syringe plunger binding or any situation requiring more force to the head of the syringe than the pump is capable of delivering. To restart the pump, remove the cause of stalling and/or increase the force setting, then press RUN. The pump will continue the interrupted procedure where it stopped.

4. Out of Range

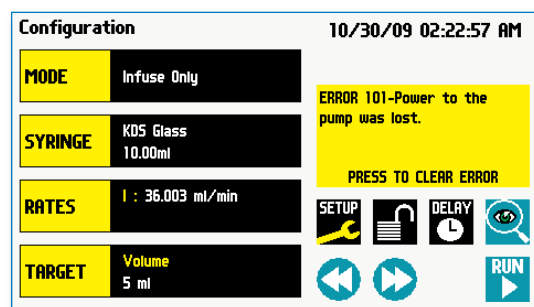
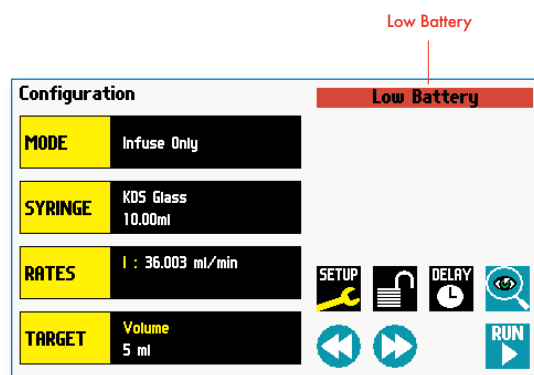
A value was entered or encountered in a pump program that was beyond the pump's limits.

5. Power Failure Notification

If power to the pump is interrupted during operation, an alarm will sound to alert the user to the interruption. When power is restored, an onscreen message is displayed to indicate the interruption.

6. Low Battery Indication

Time/Date area will turn red when battery is low.
[see Appendix E for instructions on changing the battery.]



The Cole-Parmer® Touch Screen Series Pump will notify you of a power failure with a message when power is restored

APPENDIX F: 'AUTO FILL' VALVES

There are two 'Auto Fill' Valve models available:

Item No.	Description
74905-70	Standard Pressure: Supplied with 1/8 inch ID, 1/4 inch OD silastic tubing and connector to syringe
74905-71	High Pressure: Supplied with 1/4 inch OD 304 stainless steel tubing with Swagelok fittings
74905-81	Adapter, Digital I/O

Setting Up the Valve

Set up consists of three phases: Power Connection, Pump to valve electrical connection and Valve tubing to syringe hook up.

1. Power Connection

Connect the line cord to the power supply (included). Attach the power supply to the input jack on the valve box.

2. Pump to Valve Electrical Connection

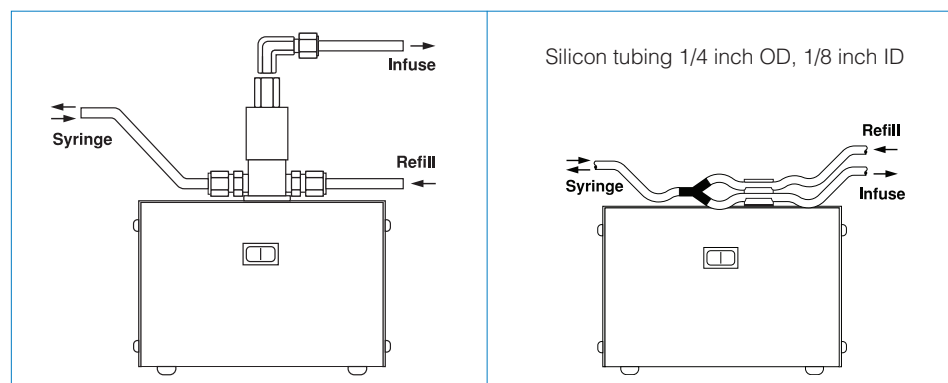
The cable fitted with a male 9 pin D-sub connector must be connected to the female D-sub connector on the rear of the Cole-Parmer® Touch Screen

Series Syringe

Pump using the Digital I/O Adapter (Order No. 74905-81). Power can now be switched on with the illuminated power On/Off switch on the front valve housing.

3. Pump to Valve Tubing Connection

Swagelok fitting should be finger tight plus 1 1/4 turn. Tubing is 304 stainless steel 1/4 inch O.D. wall thickness 0.035 inch.



APPENDIX G: ORDERING INFORMATION

Item No. Description

Model Versions

74905-50 Cole-Parmer® Touch Screen Infuse Only

74905-52 Cole-Parmer® Touch Screen Infuse/Withdraw

74905-54 Cole-Parmer® Touch Screen Infuse/Withdraw Programmable

74905-56 Cole-Parmer® Touch Screen Push-Pull

74905-58 Cole-Parmer® Touch Screen Push-Pull Programmable

Accessories

74905-60 Infuse/Withdraw 6/10 Multirack. (Six 30 to 60 ml syringes or ten 0.5 µl to 20 ml syringes) (can be sold for Infuse Only as well)

74905-61 Infuse/Withdraw 4 x 140 Multirack (Four 60 to 140 ml syringes) (can be sold for Infuse Only as well)

74905-62 Infuse/Withdraw Microliter Rack (Four 0.5 µl to 10 ml syringes) (can be sold for Infuse Only as well)

74905-70 Auto Fill Valve Box Low Pressure

74905-71 Auto Fill Valve Box High Pressure

74905-72 Anti Siphon Kit

74905-73 RS485 Pump to Pump Communication, .5 m

74905-74 RS485 Pump to Pump Communication, 2 m

74905-75 USB Cable PC to Pump Communication, 2 m

74905-76 USB Cable PC to Pump Communication, 5 m

74905-77 RS232 Cable (9 pin d-sub), 2m

74905-78 Line Cord US 115 V

74905-79 Line Cord European

74905-80 Line Cord UK

74905-81 Adapter Digital I/O (15 pin to 9 pin)

74905-82 Adapter D Sub 15 to Terminal Block

74905-83 Footswitch with Phono Jack Plug

74905-84 Adapter for 25 ml, 50 ml, 100 ml Hamilton Gastight Syringe

74905-85 Hex Key

74905-86 Lubricant SuperLube, 1 cc

Item No. Description

74905-87 Upgrade Infuse Only to Infuse/Withdraw

74905-88 Upgrade Infuse Only to Programmable

74905-89 Upgrade Infuse/Withdraw to Programmable

Option

74905-91 Internal Fan Option

74905-93 Analog Control Input Option (0 to 10 VDC)*

74905-95 Analog Control Connector

74905-96 Analog Control Cable

74905-97 Replacement Fuse

74905-98 Replacement Battery

* Only available with 74905-54, 74905-58